Reg # 211-40

Form Approved. OMB No. 2070-0060. Approval expires 2-28-95

Please read instructions on re	verse before completir	ng form.	Form Ap	proved	. OMB No. 2070-00	60. Approval expires 2-28-95		
SEPA Environmental Protection Age Washington, DC 20460					Registration	OPP Identifier Number		
			· ,		Amendment Other	220126		
Application for Pesticide - Section I								
1. Company/Product Number			2. EPA Product Mar	nager	3. F	Proposed Classification		
211-40			Mr. Marion Jo	Mr. Marion Johnson, Jr. X None Restricted				
4. Company/Product (Name) Q4.5-5.0PB-4.5			PM# 31					
5. Name and Address of Applicant (Include ZIP Code)			6. Expedited Re	6. Expedited Reveiw. In accordance with FIFRA Section 3(c)(3)				
Central Solutions, Inc. P.O. Box 15276			F	(b)(i), my product is similar or identical in composition and labeling to:				
3130 Brinke		EPA Reg. No						
Kansas C Check if this								
Check if this		Product Name						
		Se	ection - II					
Amendment - Explain	Lange Court of the	Final printed labels in repsons NOTIFICATION Agency letter dated						
Resubmission in respo	onse to Agency letter d	lated	"Me Too"	Applica	ation. ΔPR	2 4 1996		
X Notification - Explain I	oelow.		Other - Ex	plain be				
Explanation: Use addition	al page(s) if necessary	. (For section I and	Section II.) Notifica	ation	of Removal o	of Redundant Labelino		
per PR Notice 95-	2. This notif	ication is \dot{c}	onsistent w/ the	e pro	ovisions of PI	R Notice 95-2 and		
& EPA regs at 40	CFR 152.46, & 1	no other char	nges have been i	made	to the label:	ing or the CSF of		
the product. I u	nderstand that	it's a viola	ation of 18 U.S	.C. 8	Sec. 1001 to wi	illfully make any		
false statement to	o EPA. I furti 46 this produ	ner understai uct may be i	nd that if this n violation of '	not. Natra	LI. ISN'T CONS A & T may be «	subject to enforcement		
action & penaltie	s under section	ns 12 & 14 Se	ction - III of FI	FRA.	A & I likty be a	subject to enforcement		
1. Material This Product Will				- A tick over more and				
Child-Resistant Packaging	Unit Packaging	Wa	ter Soluble Packaging		2. Type of Contain	er		
X Yes	Yes		Yes		Metal	S		
No	X No	عاا	J N∘		X Plastic			
* Certification must be submitted	If "Yes" Unit Packaging wgt.		Yes" No. per kage wgt contain		Paper Other	(Specify)		
3. Location of Net Contents I	nformation	4. Size(s) Retail Cor	ntainer		cation of Label Direc	tions		
X Label C		1,2,4,8,16,3 5,6,15,30,55	2,64 oz.; 1,2.5 gallon bulk	<u>'1 </u>	on label			
6. Manner in Which Label is		X Lithograph Paper glued Stenciled		er si	lk-screened: :	s <u>leeve labe</u> ls		
			ction - IV	_		1 3		
1. Contact Point Complete	items directly below fo	or identification of in	dividual to be contacted	t, if nec	essary, to process ti	nis applicațion.)		
Name		Title			Teleph	one Not (Include Area Code)		
Paul Nobrega	and the state of t	Reg	istration Speci	alis	i	-621-6542		
		Certification				6. Data Application Received		
	y knowlinglly false or r		achments thereto are tro t may be punishable by			: \(\(\mathcal{S}\) (Stamped)		
2./Signature		3. Title	9					
Morney		No	Gistration Spe	. ما.رم	0			
4. Typed Name		5. Dat		un	120			
Paul Nobrega			March 1996					
- care womenya					<u> </u>			

Q4.5-5.0PB-4.5

DISINFECTANT - DEODORIZER - FUNGICIDE - VIRUCIDE*
WHIRLPOOL CLEANER DESIGNED FOR HOSPITAL,
NURSING HOME, AND INSTITUTIONAL USE

Active Ingredients:

n-Alkyl (60%C ₁₄ , 30%C ₁₆ , 5%C ₁₂ , 5%C ₁₈) dimethyl benzyl ammonium	chlorides	2.25%
n-Alkyl (68%C ₁₂ , 32%C ₁₄) dimethyl ethylbenzyl ammonium chlorides		2.25%
Inert Ingredients:		95.50%
-		100.00%

DANGER SEE LEFT PANEL FOR ADDITIONAL

PRECAUTIONARY STATEMENTS



P.O. BOX 15276 3130 BRINKERHOFF ROAD KANSAS CITY, KS 66115

EPA REG. NO. 211-40 EPA EST. NO. 211-KS-1

NET CONTENTS

GAI

(Left Panel)

PRECAUTIONARY STATEMENTS

Hazards to Humans and Domestic Animals

DANGER

Keep Out of Reach of Children. Corrosive. Causes severe eye and skin damage. Do not get in eyes, on skin, or on clothing. Harmful if swallowed. Wear goggles or face shield and rubber gloves when handling. Avoid contamination of food. Remove and wash contaminated clothing.

STATEMENT OF PRACTICAL TREATMENT

In case of contact, immediately flush eyes or skin with plenty of water for at least 15 minutes. For eyes, call a physician. Remove and wash all contaminated clothing before reuse. If swallowed, drink large quantities of water. Call a physician.

NOTE TO PHYSICIAN

Probable mucosal damage may contraindicate the use of gastric lavage. Measures against circulatory shock, respiratory depression, and convulsion may be needed.

ENVIRONMENTAL HAZARDS

Do not discharge effluent containing this product into lakes, streams, ponds, estuaries, oceans, or other waters unless in accordance with the requirements of a National Pollutant Discharge Elimination System (NPDES) permit and the permitting authority has been notified in writing prior to discharge. Do not discharge effluent containing this product to sewer systems without previously notifying the local sewage treatment plant authority. For guidance contact your State Water Board or Regional Office of the EPA.

STORAGE AND DISPOSAL

Do not contaminate water, food, or feed by storage or disposal.

STORAGE - Store in securely closed original container. Store in area that will prevent cross contamination of other pesticides, fertilizer, food, and feed. Avoid storage at temperature extremes or in sunlight.

PESTICIDE DISPOSAL - Pesticide wastes are acutely hazardous. Improper disposal of excess pesticides is in violation of Federal Law. If these wastes cannot be disposed of by use according to label instructions, contact your Pesticide or Environmental Control Agency, or the hazardous waste representative at the nearest EPA Regional Office for guidance.

CONTAINER DISPOSAL - Triple rinse (or equivalent), then offer for recycling or reconditioning, or puncture and dispose of in a sanitary landfill, or incineration, or if allowed by state and local authorities, by burning. If burned, stay out of smoke.

This product is not to be used as a terminal sterilant/hight level disinfectant on any surface or instrument that (1) is introduced directly into the human body, either into or in contact with the bloodstream or normally sterile areas of the body, or (2) contact intact mucous membranes but which does not ordinarily penetrate the blood barrier or otherwise enter normally sterile areas of the body. This product may be used to preclean or decontaminate critical or semi-critical medical devices prior to sterilization or high-level disinfection.

PRODUCT COMPOSITION


Q4.5-5.0PB-4.5 is a disinfectant, deodorizer, fungicide, and virucide* cleaner designed for hospital, nursing home, and institutional use. Q4.5-5.0PB-4.5 contains two highly active twin chain quaternary ammonium compounds (listed as active ingredients) for bactericidal efficacy. Inert ingredients include a balance of alkaline builders, non-ionic surfactants, and odor counteractants to enhance the cleaning and odor control characteristics. This formulation is designed to clean and disinfect bathing and therapy equipment. When used as directed, Q4.5-5.0PB-4.5 has been found effective for use on all fiberglass, stainless steel, and chrome plated surfaces of bath, therapy, whirlpools, and patient transfer equipment.

DIRECTIONS FOR USE

It is a violation of Federal Law to use this product in a manner inconsistent with its labeling.

DISINFECTION: To disinfect hard, non-porous, inanimate surfaces (such as fiberglass and stainless steel tubs and chair surfaces, chrome plated intakes and lifts, etc.), apply properly diluted Q4.5-5.0PB-4.5 with mop, cloth sponge, mechanical sprayer, or hand pump trigger sprayer so as to wet all surfaces thoroughly. For routine disinfection, proper dilution is 2 ounces of Q4.5-5.0PB-4.5 per gallon of water (1:64). Allow to remain wet for 10 minutes and then let air dry. For heavily soiled areas, a pre-cleaning step is required. Prepare a fresh dilution for each use.

FOR CLEANING FIBERGLASS BATH AND THERAPY EQUIPMENT: To remove body oils, dead tissue, hard water, and all other buildups or organic matter on inanimate surfaces after using the whirlpool unit, drain the water and refill with fresh water to just cover the intake valve. Add two ounces of Q4.5-5.0PB-4.5 for each gallon of water in the unit at this point. Briefly start the pump to circulate the solutions. Turn off pump. Wash down the unit sides, seat of the chairlift, and any/all related equipment with a clean swab or sponge. After the unit has been thoroughly cleaned, drain solutions from the unit and rinse any/all cleaned surfaces with fresh water. The unit is ready for re-use. Product to surface contact time must be at least 10 minutes for proper disinfection.

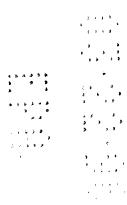



PERFORMANCE DATA

Efficacy tests have demonstrated that Q4.5-5.0PB-4.5 is an effective bactericide, fungicide, and virucide* on hard, non-porous, inanimate surfaces in the presence of moderate amounts of organic soil (5% blood serum). All testing has been performed using A.O.A.C. protocol testing as required by the EPA at a use-dilution of 1:64 (two ounces per gallon of water) in the presence of 5% blood serum (used to simulate organic soil load), at 20 degrees celcius with an exposure time of 10 minutes.

GERMICIDAL

The broad spectrum activity of Q4.5-5.0PB-4.5 has been confirmed against the following organisms:

Pseudomonas aeruginosa

Escherichia coli

Serratia marscenscens

Proteus vulgaris

Methicillin resistant Staphylococcus aureus (MRSA)

Staphylococcus aureus Klebsiella pneumoniae Salmonella choleraesuis Candida albicans

FUNGICIDAL

Q4.5-5.0PB-4.5's efficacy as a fungicide has been confirmed through testing against the pathogenic fungi, Trichophyton mentagrophytes.

VIRUCIDAL

*Testing of Q4.5-5.0PB-4.5 has confirmed its virucidal activity against many lipophilic viruses including:

Influenza A2-Asian......(representative of the common flu virus)

Herpes simplex Type 1...(causative agent of fever blisters)

Adeno virus type 2......(causative agent of upper respiratory infections)

Vaccinia virus.....(representative of the pox virus group)

KILLS HIV-1 ON PRECLEANED ENVIRONMENTAL SURFACES/OBJECTS PREVIOUSLY SOILED WITH BLOOD/BODY FLUIDS in health care settings or other settings in which there is an expected likelihood of inanimate surfaces/objects likely to be soiled with blood or body fluids, and in which the surfaces/objects likely to be soiled with blood or body fluids can be associated with the potential for transmission of human immunodeficiency virus Type-1 (HIV-1) (associated with AIDS). Contact time: ten (10) minutes.

SPECIAL INSTRUCTIONS FOR CLEANING AND DECONTAMINATION AGAINST AIDS (HIV-1 VIRUS) ON SURFACES/OBJECTS SOILED WITH BLOOD/BODY FLUIDS WITH TEN (10) MINUTE CONTACT TIME.

PERSONAL PROTECTION: Disposable gloves should be worn when cleaning blood, body fluids, or contaminated items. Hands should be washed immediately and thoroughly after gloves are removed, or if hands are contaminated with blood or body fluids. Masks and protective eyewear or faceshields should be worn during cleaning procedures involving blood or body fluids. Gowns or aprons should be worn during leaning procedures for protection against splashes of blood or other body fluids.

CLEANING PROCEDURE: Blood and other fluids must be thoroughly cleaned from surfaces and objects before application of the disinfectant. After cleaning is complete, follow disinfection described previously.

DISPOSAL OF INFECTIOUS MATERIALS: Contaminated waste and blood and other body fluids should be autoclaved and disposed according to local regulations for infectious waste disposal.