

Front

PU 22
1/5-1/2
11/16

**BULK
LYSOL®
BRAND DISINFECTANT**

GERMICIDAL* Staphylocidal, Pseudomonacidal**, Tuberculocidal
FUNGICIDAL* (Against pathogenic fungi)
VIRUCIDAL* (Against Influenza Virus Type A₂ [Japan], Herpes Simplex Type 1, Vaccinia, Adenovirus Type 2 and Human Immunodeficiency Virus Type 1 (HIV-1) (AIDS Virus))

*On Hard, Nonporous Surfaces **At 2% Solution

Prevents Mold and Mildew

ACTIVE INGREDIENTS: o-Phenylphenol, 2.80%, o-Benzyl-p-chlorophenol, 2.70%.

INERT INGREDIENTS: 94.50%.

EPA REGISTRATION NO. 675-19 EPA EST. NO. 777-OH-1

KEEP OUT OF REACH OF CHILDREN

WARNING: Harmful if swallowed. Causes eye and skin irritation. See Back Panel for Statement of Practical Treatment and Additional Precautionary Statements.


NET CONTENTS ONE GALLON

Alternate Sizes

**NET CONTENTS: 64 FL. OZ. (1.892 liters)
5 GALLONS
55 GALLONS**

5.4.90

Not Reviewed. Registration to be in accordance with
the Federal Insecticide, Fungicide, and Rodenticide Act
(FIFRA) and the Federal Food, Drug, and Cosmetic Act
(FDCA).

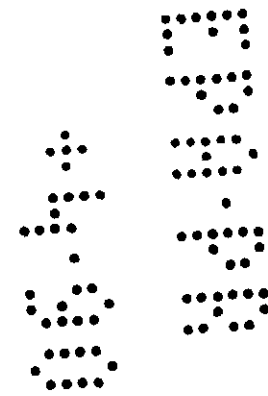

Pg 200

HOW TO USE BULK LYSOL® BRAND DISINFECTANT
(See Brochure for detailed instructions)

Bulk LYSOL Brand Disinfectant is recommended for use in hospitals, nursing homes, dental offices and other institutional facilities. Bulk LYSOL Brand Disinfectant may be used on floors, walls, bathroom fixtures and surgical instruments. (See brochure for additional information).

BULK LYSOL IS IN CONCENTRATED FORM. DILUTE ACCORDING TO DIRECTIONS. It is a violation of federal law to use this product in a manner inconsistent with its labeling.

To prepare proper dilutions in required quantities, add full strength Bulk LYSOL to water as below.

FOR CLEANING OF SURGICAL AND DENTAL INSTRUMENTS:

Prior to disinfection instruments must be cleaned to remove all adhering blood and serous exudates, body fluids and human wastes.

FOR DISINFECTION OF SURGICAL AND DENTAL INSTRUMENTS:


Immerse instruments in Bulk LYSOL Brand Disinfectant solution for 15 minutes to wet all surfaces thoroughly. Flush and fill the lumen of hollow instruments/equipment, avoiding entrapment of air bubbles. If instruments are to remain in solution for substantially longer periods of time, add 1 gram of sodium bicarbonate to each quart of solution; stir to dissolve completely. Also add 2 oz. of Isopropyl alcohol per quart to solution. Fresh solution should be prepared at least daily or more often if the solution becomes diluted or soiled.

Correct
Solution
Strength

BULK
LYSOL WATER

3%
(3:100)

30 cc.	1 quart
4 oz.	1 gallon
12 oz.	3 gallons
20 oz.	5 gallons

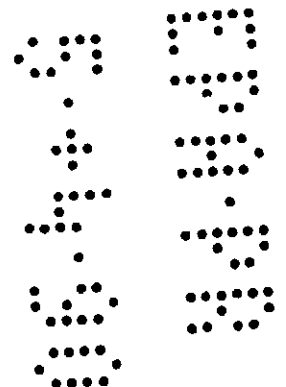

02.3016


**FOR SANITARY DISPOSAL OF
URINE, FECES:**

To urine, add equal amount of Bulk LYSOL Brand Disinfectant solution. For feces, place solution in pan in amount slightly larger than expected discharge. After receiving feces, mix thoroughly with solution then cover and let stand for one hour before disposal.

**FOR SANITIZING LAUNDRY, LINENS,
BLANKETS AND DIAPERS:**

Add 1 cup of Bulk LYSOL Brand Disinfectant to each 17 gallons of water.

(See brochure for detailed instructions)

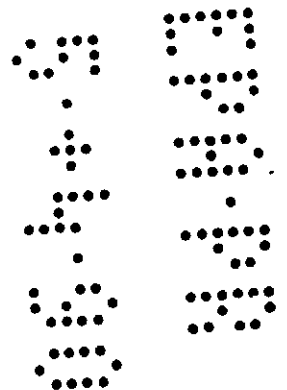

**BULK
LYSOL®
BRAND DISINFECTANT**

PRECAUTIONARY STATEMENTS: HAZARDS TO HUMANS AND DOMESTIC ANIMALS

WARNING: Harmful if swallowed. Causes eye and skin irritation. Do not get in eyes. Avoid contact with eyes and skin. Avoid contamination of food.

STATEMENT OF PRACTICAL TREATMENT:

In case of eye contact, **IMMEDIATELY** flush eyes thoroughly with water, remove any contact lenses, and continue to flush eyes with plenty of water for at least 15 minutes. Get medical attention.

In case of skin contact, wash with soap and water. If irritation persists, consult a physician.


If swallowed, drink promptly a large quantity of water to dilute the product. Get immediate medical attention.

STORAGE AND DISPOSAL

64 fl. oz.
& 1 gallon **STORAGE:** Store in original container in areas inaccessible to small children. Keep securely closed.
 DISPOSAL: Do not reuse empty container. Rinse thoroughly with water and detergent. Wrap and discard in trash.

5 & 55
gallons Do not contaminate water, food, or feed by storage and disposal.
 STORAGE: Store in original container in areas inaccessible to small children. Keep securely closed.
 Pesticide Disposal: Wastes resulting from the use of this product may be disposed of on site or at an approved waste disposal facility.
 Container Disposal: Triple rinse (or equivalent). Then offer for recycling or reconditioning, or puncture and dispose of in a sanitary landfill, or incineration, or, if allowed by state and local authorities, by burning. If burned, stay out of smoke.

**NATIONAL LABORATORIES
Lehn & Fink Products Group
Montvale, NJ 07645**

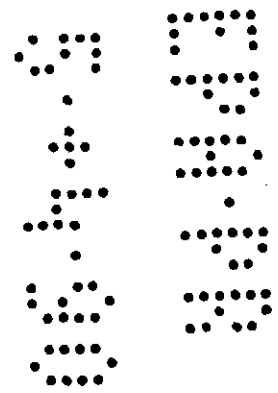

HOW TO USE BULK LYSOL® BRAND DISINFECTANT


(See Brochure for detailed instructions)

BULK LYSOL IS IN CONCENTRATED FORM. DILUTE ACCORDING TO DIRECTIONS.

	Correct Solution Strength	To prepare proper dilutions in required quantities, add full strength Bulk LYSOL to water as below.	
		BULK LYSOL	WATER
FOR PREVENTING RUST IN HEAT DISINFECTION:			
Bulk Lysol Brand	1/2% (1:200)	20 cc. 2 oz.	1 gallon 3 gallons
Disinfectant neutralizes the rusting action of boiling water.			

FOR GENERAL DISINFECTION:(including areas occupied by
tubercular patients)To disinfect hard non-porous
surfaces, first clean surface
thoroughly prior to applica-
tion of Bulk LYSOL Brand
Disinfectant. Apply solution
to surface using a sponge or
mop to wet all surfaces
thoroughly. Letstand 10 minutes before
wiping or allow to air dry. 1%
(1:100)Fresh solution should be
prepared daily or more often
if solution becomes diluted or
soiled. For tuberculocidal
activity 1% Bulk LYSOL Brand
Disinfectant is effective
in 10 minutes at 20°C (68°F).**FOR MOLD AND MILDEW ON HARD
NON-POROUS SURFACES:**Wash surfaces with Bulk LYSOL
solution; allow to air dry.
Repeat application as
necessary.

10 cc.	1 quart
40 cc.	1 gallon
4 oz.	3 gallons

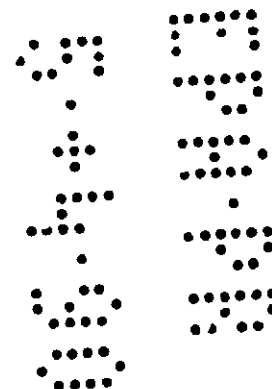FOR PSEUDOMONACIDAL ACTIVITY:Clean surface thoroughly. Apply
solution using a sponge or mop.
Wet surface thoroughly. Let
stand 10 minutes before wiping
or let it air dry.

Panel 4 (continued)

PURPOSE

FOR HIV-1 (AIDS VIRUS) ACTIVITY:
KILLS HUMAN IMMUNODEFICIENCY
VIRUS TYPE 1 (HIV-1) (AIDS
VIRUS) ON PRECLEANED
ENVIRONMENTAL SURFACES/OBJECTS
PREVIOUSLY SOILED WITH BLOOD/
BODY FLUIDS.

SPECIAL INSTRUCTIONS FOR CLEANING
AND DECONTAMINATION AGAINST HIV-1
OF SURFACES/OBJECTS SOILED WITH
BLOOD/BODY FLUIDS.

PERSONAL PROTECTION: When
handling items soiled with blood
or body fluids, use disposable
latex gloves, gowns, masks or
eye coverings.

CLEANING PROCEDURES: Blood and
other body fluids must be
thoroughly cleaned from surfaces
and objects before application of
Bulk LYSOL Brand Disinfectant.

CONTACT TIME: Allow surface
to remain wet for 10 minutes.

DISPOSAL OF INFECTIOUS MATERIALS:
Blood and other body fluids should
be autoclaved and disposed of
according to local regulations for
infectious waste disposal.

Correct
Solution
Strength
For HIV-1
Activity

2 $\frac{1}{2}$
(2:100)

To prepare proper dilutions
in required quantities,
add full strength
Bulk LYSOL to water
as below.

BULK

LYSOL

80 cc.

8 oz.

WATER

1 gallon

3 gallons

ONE TABLESPOONFUL EQUALS ABOUT 1/2 FLUID OUNCE OR 15 CC.

