

PU 32.
6/1/54
ME

FRONT PANEL

O-SYL
DISINFECTANT . DETERGENT

GERMICIDAL * Staphylocidal * Pseudomonacidal *Tuberculocidal
FUNGICIDAL * (Against pathogenic fungi)
VIRUCIDAL * (Against Influenza A₂, Herpes simplex, Vaccinia, and
Adenovirus Type 2)
*On Environmental Surfaces

ACTIVE INGREDIENTS:

o-Phenylphenol5.00%
o-Benzyl-p-chlorophenol4.50%

INERT INGREDIENTS:90.50%
TOTAL..100.00%

EPA Registration. No. 675-24

DANGER: KEEP OUT OF REACH OF CHILDREN
Corrosive, Causes eye damage and skin irritation. Do not get in
eyes, on skin or on clothing. Protect eyes and skin when handling
concentrate. Harmful or fatal if swallowed. Avoid contamination
of food. Rinse empty container thoroughly with water and discard
it.


See First Aid statement and other precautions on side panel.

NET CONTENTS: ONE GALLON

National Laboratories
Lehn & Fink Products Group
Montvale, NJ 07645 Plant: Toledo, Ohio 43612

5,460

Not Reviewed. Registrant claims to be in accordance with
Paragraph 101 of Part 161 of the Code of Federal Regulations for Phase 2
Registration. Efficacy of this product is subject to change in Status.

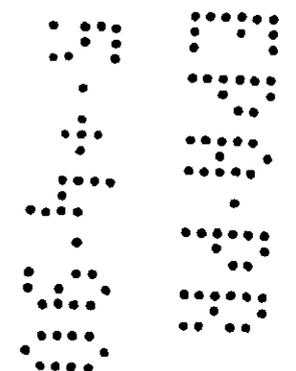

Right Side Panel

HOW TO USE HOSPITAL TYPE
O-SYL
DISINFECTANT . DETERGENT

O-SYL IS IN CONCENTRATED FORM, DILUTE ACCORDING TO DIRECTIONS.
ONE TABLESPOON EQUALS ABOUT 1/2 FLUID OUNCE OR 15 CC.

To prepare proper dilutions in required quantities, add full strength O-SYL to water as below:


PURPOSE	Correct Solution Strength	O-SYL	WATER
FOR GERMICIDAL CLEANING AND DEODORIZING. (Including premises occupied by tubercular patients) Floors, Walls, Furniture, Food and linen cars, Garbage cans.		8 cc.	1 quart
		30 cc.	1 gallon
	0.78%	3 oz.	3 gallons
	(1:128)	5 oz.	5 gallons
FOR DISINFECTION OF: Thermometers, surgical and dental instruments.	3 %	30 cc.	1 quart
	(3:100)	4 oz.	1 gallon
		12 oz.	3 gallons
		20 oz.	5 gallons

See brochure for detailed instructions.

FIRST AID:

In case of contact, immediately flush eyes or skin with plenty of water for at least 15 minutes. For eyes, call a physician. Remove and wash contaminated clothing before reuse.

If swallowed, drink promptly a large quantity of milk, egg whites, or gelatin solution. If these are available drink large quantities of water. Avoid alcohol. Call a physician immediately.

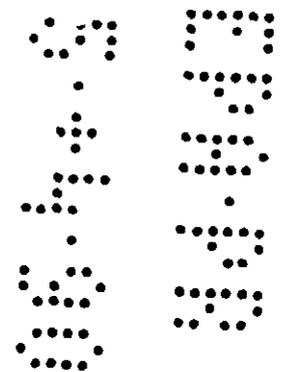

Left Side Panel

AMERICA'S FOREMOST MANUFACTURER OF
DISINFECTANTS AND DISINFECTANT-DETERGENTS
FOR ENVIRONMENTAL SANITATION

