

8-108-81

Burnishine Products  
8140 North Ridgeway Avenue  
Skokie, IL 60076

Attention: Carl B. DeMasi

Gentlemen:

Subject: Burnishine Germicidal Cloth  
EPA Registration No. 1130-6  
Your Amendment Application Dated September 10, 1986

The amendment referred to above, submitted in connection with registration under section 3(c)(7)(A) of the Federal Insecticide, Fungicide, and Rodenticide Act (FIFRA), is acceptable provided that you:

1. Submit/cite all data required for registration/reregistration of your product under FIFRA section 3(c)(5) when the Agency requires all registrants of similar products to submit such data.
2. Make the labeling changes listed below before you release the product for shipment bearing the amended labeling:
  - a. Packaging this product as an individual towelette and as multiple packs in a cylindrical reclosable tub would be acceptable.
  - b. However, for the multiple packs in addition to the "Net Wt" the label should give the "Net Contents" and state the number of towelettes.
3. Submit five (5) copies of your final printed labeling before you release the product for shipment.

90163:Lee:L-5:KENCO:12/10/86:12/22/86:eg:VO

CONCURRENCES								
SYMBOL ▶								
SURNAME ▶								
DATE ▶								


2/8

2

If these conditions are not complied with, the registration will be subject to cancellation in accordance with FIFRA section 6(e). Your release for shipment of the product constitutes acceptance of these conditions.

A stamped copy of the label is enclosed for your records.

Sincerely yours,

John H. Lee   
Product Manager (31)  
Disinfectants Branch  
Registration Division (OS-767C)

Enclosure

## MASTER LABEL COPY

Front Copy

BURNISHINE<sup>R</sup> GERMICIDAL CLOTH

For Use In: CRITICAL CARE AREAS

Where control of the hazards of cross-contamination is of prime importance.

Effective against a wide range of bacteria including Staphylococcus aureus, Salmonella choleraesuis and Pseudomonas aeruginosa.

A novel way to effectively disinfect environmental surfaces.

Aids in the reduction of cross-contamination.

Convenient -- Pre-Measured -- Disposable

## ACTIVE INGREDIENTS

n-Alkyl (68% C <sub>12</sub> , 32% C <sub>14</sub> )	
dimethyl ethylbenzyl ammonium chlorides .....	0.14%
n-Alkyl (60% C <sub>14</sub> , 30% C <sub>16</sub> , 5% C <sub>12</sub> , 5% C <sub>18</sub> )	
dimethyl benzyl ammonium chlorides .....	0.14%
Isopropyl alcohol .....	8.00%

INERT INGREDIENTS .....	91.72%
-------------------------	--------

Keep Out of Reach of Children

CAUTION

See additional precautions on back panel

NET WT.

EPA Reg. No. 1130-6

EPA Est. No. 10036NJ-1

JAN 08 1967

Under the Federal Insecticide, Fungicide, and Rodenticide Act, as amended, for the pesticide registered under EPA Reg. No.

1130-6

4/8 .

Back Panel

A dual chain quaternary/alcohol solution impregnated in a wiping cloth.

- \* Assures proper dilution strength and product effectiveness and standardized technique.
- \* Greater quaternary dispersal with no hazard of cross-contamination quick drying.
- \* Non-woven disposable cloth.
- \* Effective against Staphylococcus aureus, Salmonella choleraesuis, Pseudomonas aeruginosa.

For use in hospitals and other critical care areas where control of the hazards of cross-contamination is required. Use on surfaces and equipment such as tables, carts, baskets, counters, cabinets, telephones in operating rooms, intensive care and emergency units. ;

Use on stainless steel, Formica, glass and other hard, non-porous surface.

Some areas of use: Hard surfaces in: Surgery, Recovery, Anesthesia, X-Ray Cat. lab, E.R., Orthopedics, New Born Nursery, Respiratory Therapy, Radiology, Central Supply.

DIRECTIONS FOR USE: It is a violation of Federal law to use this product in a manner inconsistent with its labeling.

TOWELETTE

Remove gross filth. Open, unfold, wipe to thoroughly wet surface. Allow to remain wet 10 minutes, let air dry.

DISPENSER

To start feed: Flip open cap. Remove and discard foil seal. remove large cover from container. From center of towelette roll, pull up a towelette corner, twist it into a point and thread it through the hole in the container cover. Pull through about one inch. Replace cover. Pull out first towelette and snap off at a 90° angle. Remaining towelettes feed automatically, ready for next use. When through using, keep small center cap closed to prevent moisture loss.

Remove gross filth. Wipe to thoroughly wet surface. Allow to remain wet 10 minutes, let air dry.

JAN 08 1961

UNITED STATES  
DEPARTMENT OF  
HEALTH, EDUCATION  
AND WELFARE  
REGISTERED UNDER LHA Reg. No.

1130-6

Back Panel, Cont'd.

PRECAUTIONARY STATEMENTS Hazards to humans and domestic animals.

CAUTION Avoid contact with eyes. In case of eye contact, flush with plenty of water for at least 15 minutes. If irritation persists, call a physician.

STORAGE AND DISPOSAL

TOWELETTE

Do not reuse towel. Dispose of used towel and empty packet in trash.

DISPENSER

Do not reuse towel or empty container. Wrap container and put in trash collection. Dispose of used towels in trash.


BURNISHINE PRODUCTS

DIV. OF HERBERT STANLEY COMPANY \* SKOKIE, ILLINOIS 60076

JAN 08 1967

U. S. Environmental Protection Agency,  
Washington, D. C. 20460  
is amended, for the purpose  
registered under EPA Reg. No.

1130-6

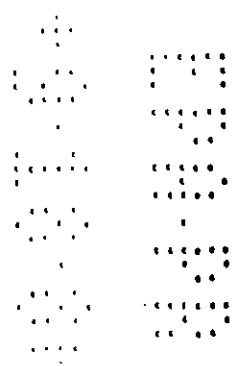

6/8 .  
Consumer Cloth  
Front Side

1 Step Cleaner  
& Disinfectant, Sanitizer (Non-Food - Contact Surfaces)  
Deodorizer for Disinfecting Contaminated Surfaces

A novel way to instantly clean, deodorize and sanitize against odor-causing organisms and disinfect hard surfaces. Aids in reducing cross-contamination. Kills dangerous germs such as Salmonella choleraesuis, Staphylococcus aureus and Pseudomonas aeruginosa. Easy to use. Disposable. For the reasonable person who wants to be sure. For use in public toilets, laundromats, sick rooms, hospitals, camp sites, RV's, motels, factories and offices -- whenever sanitation is of prime importance.

JAN 08 1967  
Under the President John F. Kennedy  
Executive Order 11652, for the purpose of  
disinfecting under LPA Reg. No.  
1130-6

7/8  
Back Side  
Consumer Cloth

This product kills a wide range of organisms including gram negative and gram positive bacteria such as *Pseudomonas aeruginosa*, *Staphylococcus aureus*, *Salmonella choleraesuis* which cause such diseases as boils, eye and skin infections and food poisoning.

For use in households; in hospitals - veterinary clinics; in laundromats in washers and dryers; on folding tables; on toilet seats and bathroom fixtures; hampers and bedpans; sickrooms; offices; on telephones; and hard hats and headphones. For use in dairies, food processing plants, food establishments, restaurants, taverns and commercial kitchens (disinfected food contact surfaces must be rinsed with potable water before reuse); in hotels, in public restroom accommodations and commercial laundries; on stainless steel, chrome, porcelain, glass, formica, vinyl, plastic and other hard, non-porous surfaces and CPR training manikins.

#### DIRECTIONS FOR USE

It is a violation of Federal Law to use this product in a manner inconsistent with its labeling.

#### To Clean, Sanitize (Non-Food Contact Surfaces) And Deodorize

TOWELETTE	Remove gross filth. Open, unfold, wipe until surface is thoroughly wet. Allow to remain wet 5 minutes. Dry with tissue or hand towel.
DISPENSER	To start feed: Flip open cap. Remove and discard foil seal. Remove large cover from container. From center of towelette roll, pull up a towelette corner, twist it into a point and thread it through the hole in the container cover. Pull through about one inch. Replace cover. Pull out first towelette and snap off at a 90° angle. Remaining towelettes feed automatically, ready for next use. When through using, keep small center cap closed to prevent moisture loss. Remove gross filth. Wipe to thoroughly wet surface. Allow to remain wet 5 minutes. Dry with tissue or hand towel.

8/4 .  
Back Panel, Cont'd.

To Disinfect

TOWELETTE

Remove gross filth. Open package, remove towel and wipe surface with towel until surface is thoroughly wet. Allow to remain wet 10 minutes, let air dry.

DISPENSER

To start feed: Flip open cap. Remove and discard foil seal. Remove large cover from container. From center of towelette roll, pull up a towelette corner, twist it into a point and thread it through the hole in the container cover. Pull through about one inch. Replace cover. Pull out first towelette and snap off at a 90° angle. Remaining towelettes feed automatically, ready for next use. When through using, keep small center cap closed to prevent moisture loss. Remove gross filth. Wipe to thoroughly wet surface. Allow to remain wet 10 minutes, let air dry.

ADDITIONAL PRECAUTIONS

Caution: Avoid contact with eyes. In case of eye contact, flush with water. If irritation persists, get medical attention.

STORAGE & DISPOSAL

TOWELETTE

Do not reuse towel. Dispose of towel and empty packet in trash.

DISPENSER

Do not reuse towel or empty container. Wrap container and put in trash collection. Dispose of used towels in trash.