

Label

Front Panel

0.15 Fl.Oz.

THERMO BENZ

1:750 Tincture Benzalkonium Chloride
with Alcohol 50%

Active Ingredients:

n-Alky1 (C ₁₂ 50%, C ₁₄ 30%, C ₁₆ 17%, C ₈₋₁₈ 3%) Dimethyl Benzyl Ammonium Chloride (Benzalkonium Chloride)	0.13%
Alcohol	50 %

Inert Ingredient	49.87%
------------------	--------

CAUTION: Keep out of the reach of children. Avoid contamination
of food.

XTTRIUM LABORATORIES, INC.
CHICAGO, ILLINOIS 60609


Back Panel

Disinfectant
for the Storage of Fever Thermometers

DIRECTIONS: Remove and discard cap. Place vial in wall clip.
Insert thermometer. Change vial with each new
patient or every seven days. To remove vial
from wall clip pull forward from tape.

USDA Reg. # _____

CONTROL _____

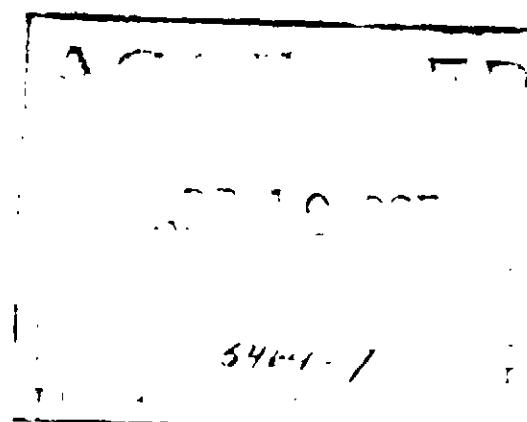

Formula

THERMO BENZ

Per Vial

Per 100 gal. Batch

Benzalkonium Chloride USP

6 mg.

1 lb. 1 oz. 351 grains

Alcohol USP

2.36 cc

52 gal. 74 oz.

Water, Purified USP

q.s.4.5 cc

q.s. 100 gal.