

# Carrier Chemicals Inc.

Main Office and Plant: VERNON, N. J. 07462

Library Copy


PAGE 1 (1)

Use CONQUEST 100 for

MAIL ADDRESS

CARRIER CHEMICALS INC.  
VERNON, NEW JERSEY 07462

All Floors  
Walls  
Painted Surfaces  
Laminate Surfaces  
Isolation Areas  
Surgical Areas  
Cosmetical Areas  
Nurseries  
Laboratories

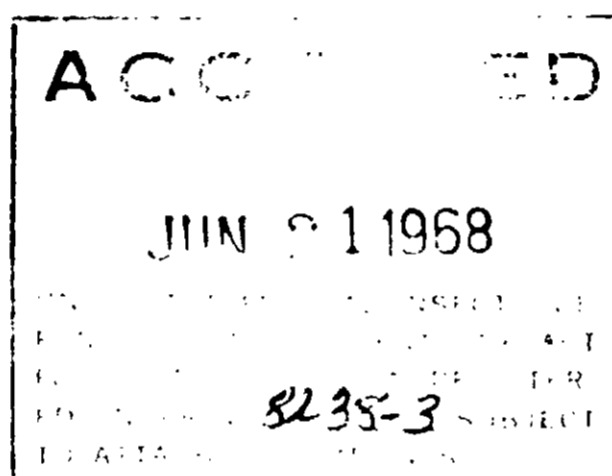

Instruments  
Utensils  
Basins  
Tubing  
Catheters  
All Plastics  
All Rubber  
Light Fixtures  
Cabinet Tops

**DIRECTIONS FOR USE:**

For the general cleaning and sanitizing of walls, floors, utensils and equipment. For floor mopping and dusting use 1 1/4 oz. Conquest 100 per gallon of water (1.25% solution). Allow the solution a few minutes of contact with surfaces being cleaned.

To achieve the maximum degree of safety in highly contaminated areas such as Isolation rooms, use 2 oz. Conquest 100 per gallon of water (1.5% solution).

The disinfection of walls and floors is not necessary.

Conquest 100 is harmless to resilient floors, painted and varnished surfaces, tile and terrazzo. This solution should never be mixed with other chemicals, its activity can be altered by such mixing.

**CAUTION:**

BE CAREFUL IF SWALLOWED. Avoid prolonged skin contact and inhalation in gases. In case of contact of concentrate with skin, wash by flooding with large amounts of water. Also skin surface with soap and water.



**Carrier** *Chemicals Inc.*

Main Office and Plant: VERNON, N. J.

NET  
GROSS  
8-17-77

PROPOSED LABEL PANEL 1

USE CONQUEST 100 FOR

All Floors	Instruments
Walls	Utensils
Painted Surfaces	Basins
Waxed Surfaces	Tubing
Isolation Areas	Catheters
Surgical Areas	All Plastics
Obsterical Areas	All Rubber
Nuseries	Light Fixtures
Bathrooms	Counter Tops
Toilets	Thermometers
Laundries	Waste Cans

DIRECTIONS FOR USE

For the general cleaning and disinfecting of walls, floors, utensils and equipment; for damp mopping and dusting use 1 1/4 oz. CONQUEST 100 per gallon of water (a 1% solution). Allow the solution a few minutes of contact with surfaces being cleaned.

To achieve the maximum margin of safety in heavily contaminated areas, and for heavy duty cleaning operations use 2 oz. CONQUEST 100 per gallon of water (1.5% solution). The rinsing of walls and floors is not necessary.

When used as directed CONQUEST 100 is harmless to fabrics, resilient floors, painted and varnished surfaces, tile and terrazzo. This formula should never be mixed with other chemicals, its activity can be altered by such mixing.

CAUTION

Harmful if swallowed. Avoid prolonged skin contact and splashing in eyes. In case of accidental contact of concentrate with skin wash by flooding with large amounts of water. Wash skin surfaces with soap and water.



PROPOSED LABEL PANEL 2

PHENOLIC DETERGENT DISINFECTANT

CONQUEST 100

DEODORIZES AND DISINFECTS WHILE IT CLEANS- IN ONE OPERATION

Broad spectrum germicidal concentrate and deodorant with an effective detergent system, designed for use in hot or cold, and hard or soft water.

Phenol Coefficient ( A.O.A.C. )

S. aureus - 9  
S. typhosa - 9

Use-Dilution Conformation ( A.O.A.C. )

S. aureus - 1:100  
S. choleraesuis - 1:100  
A. aerogenes - 1:100  
E. coli - 1:100  
T. interdigitalis - 1:100  
S. fecalis - 1:100  
P. aeruginosa - 1:90  
S. typhosa - 1:100

ACTIVE INGREDIENTS:

Isopropyl Alcohol	10.00 %
Sodium 4-chloro-2-cyclohexyl phenol	3.00 %
Sodium ortho-benzyl para-chloro phenol	0.75 %
Tetrasodium ethylene diamine tetracetate	0.50 %
	<hr/>
	14.25 %
INERT INGREDIENTS:*	85.75 %

\*Active as cleaning agents.

CAUTION

KEEP OUT OF REACH OF CHILDREN

U.S.P.A. Registration No. \_\_\_\_\_

CARRIER CHEMICALS, INCORPORATED

Main Office and Plant: Vernon, New Jersey 07462



# Carrier Chemicals Inc.

Main Office and Plant: VERNON, N. J.

MAIL ADDRESS  
 BOX 380 R. D. 1  
 WARWICK NEW YORK

PROPOSED LABEL PANEL 3

DILUTION TABLE

<u>SOLUTION STRENGTH</u>	Quantity of CONQUEST 100 to make:				
	<u>1 qt.</u>	<u>1 gal.</u>	<u>2 gal.</u>	<u>3 gal.</u>	<u>5 gal.</u>
1.0%	10 cc.	1 1/4 oz.	2 1/2 oz.	4 oz.	6 1/4 oz.
1/5%	1/2 oz.	2 oz.	4 oz.	6 oz.	10 oz.
2.0%	20 cc.	2 1/2 oz.	5 oz.	8 oz.	12 1/2 oz.
3.0%	1 oz.	4 oz.	8 oz.	12 oz.	20 oz.

SPECIAL APPLICATIONS

INSTRUMENTS: Heat sensitive instruments can be cold disinfected by soaking in a 3% solution of CONQUEST 100 (1 oz. per quart of water). Follow by a clear water rinse.

ROEPPER and PLASTIC GOODS: For the disinfection of catheters, gloves, masks, tubes and other items by the cold disinfection method soak for exactly 15 minutes (IMPORTANT!!) in a 3% CONQUEST 100 solution. Rinse with water. Catheters should also be soaked for one minute in 70% alcohol as a final step.

SPECIALIZED CLEANING TECHNIQUES: Consult your Carrier representative for information regarding other specialized cleaning and disinfecting operations.