

8383-7

05/13/2010

1/4

UNITED STATES ENVIRONMENTAL PROTECTION AGENCY



UNITED STATES ENVIRONMENTAL PROTECTION AGENCY
WASHINGTON, DC 20460

MAY 13 2010

OFFICE OF
PREVENTION, PESTICIDES
AND TOXIC SUBSTANCES

Betsy Potter
CONTEC, Inc.
525 Locust Grove
Spartanburg, SC 29303

Subject: **Sporicidin Disinfectant Towelettes**
EPA Registration No.: 8383-7
Application Date: May 13, 2010
Receipt Date: May 13, 2010

Dear Ms. Potter:

This acknowledges receipt of your notification, submitted under the provision of PR Notice 98-10, FIFRA section 3(c) 9.

Proposed Notification:

- Adding Text: Fresh Scent descriptor
- Modify Confidential Statement of Formula (CSF) to add cleared fragrances

General Comment:

Based on a review of the material submitted, the following comment applies:

A copy of the accepted notification and CSF has been inserted in your file for future reference.

Should you have any questions or comments concerning this letter, you may contact me by telephone at (703) 308-6416 or by e-mail at campbell-mcfarlane.Jacqueline@epa.gov or Stacey Grigsby by telephone at (703) 305-6440 or by e-mail at grigsby.stacey@epa.gov during the hours of 8:00am to 4:00pm EST. When submitting information or data in response to this letter, a copy of this letter should accompany the submission to facilitate processing.

Sincerely,

Jacqueline McFarlane
(Acting) Product Manager (34)

CONCURRENTLY Antimicrobials Division (7510P)

SYMBOL							
SURNAME							
DATE							



March 22, 2010

Ms. ShaRon Carlisle, Acting PM # 34
Document Processing Desk (NOTIF)
Office of Pesticide Programs (7504P)
U.S. Environmental Protection Agency
1200 Pennsylvania Ave., NW
Washington, DC 20460-0001

**Subject: Contec, Inc.
Sporicidin Disinfectant Towelettes, EPA Reg. 8383-7, Notification of label change per PR Notice PR
Notice 98-10**

Dear Ms. Carlisle:

This is a Notification of label change per PR Notice 98-10 to add the text "Fresh Scent" to the label for 8383-7, in accordance with the recently accepted alternate fragrance formulation for this product.


Enclosed please find form 8570-1 and a copy of the label with changes marked in red.

This notification is consistent with the guidance in PR Notice 98-10 and the requirements of EPA's regulations at 40 CFR §§ 156.10, 156.140, 156.144, 156.146, and 156.156. No other changes have been made to the labeling or the Confidential Statement of Formula for this product. I understand that it is a violation of 18 U.S.C. Sec. 1001 to willfully make any false statement to EPA. I further understand that if the amended label is not consistent with the requirements of 40 CFR §§ 156.10, 156.140, 156.144, 156.146, and 156.156, this product may be in violation of FIFRA and I may be subject to enforcement action and penalties under sections 12 and 14 of FIFRA.

Thank you for your assistance with this submission. If you have any questions, please call me at 864-699-8321 or e-mail me at bpotter@contecinc.com.

Sincerely,

Betsy Potter
Contec, Inc.

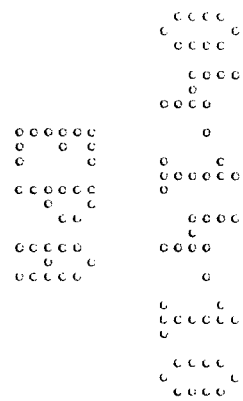



CONTINUOUS RESIDUAL ACTIVITY
UP TO 6 MONTHS*

Sporicidin[®]

BRAND

DISINFECTANT TOWELETTES

An effective disinfectant and active deodorant for hospitals, clinics, medical, dental and veterinary offices, farm premises, laboratories, industrial clean rooms, homes, nursing homes, ambulances, hotels, restaurants, schools, airplanes, trains, boats, autos, buses, health spas and toilets.

Active Ingredients
Phenol1.56%
Sodium Phenate0.06%
Inert Ingredients98.38%
.....100.00%

Out of Reach of Children

CAUTION

Net contents: 180 Towelettes 5" x 8" saturated with disinfectant; 1 lb. 11 oz.

- BACTERICIDAL
- FUNGICIDAL
- VIRUCIDAL**
- TUBERCULOCIDAL
- PREVENTS MOLD AND MILDEW
- ELIMINATES ODORS*

EPA REG. NO. 8383-7
EPA Est. Nos. 034346-NB-001,
08255-IL-001, 71979-SC-001

See batch code for actual Est. No.

Part No.: CAN-18012F



INSTITUTIONAL • FOOD PROCESSING • RESIDENTIAL

AN EFFECTIVE CLEANER, DISINFECTANT AND ACTIVE DEODORIZER FOR INSTRUMENTS, AND EQUIPMENT ON PRE-CLEANED, HARD, NON-POROUS surfaces in health/hospital treatment and patient rooms, operating rooms, schools, ambulances, hotels, restaurants, boats, planes, buses, and the home. Recommended for medical and dental equipment, beds, surgical carts, counter tops, mannequins, telephones, keyboards, furniture, wheelchairs, walkers, sinks, floors, walls, light switches, linen hampers, shower stalls, toilet seats, trash containers and damp and musty areas, including animal areas. Self-sanitizes against odor-causing bacteria and fungi. Reduces the hazards of cross-contamination.

PASSES AOAC EFFICACY STANDARDS FOR DISINFECTION IN THE PRESENCE OF 5% ORGANIC SOIL. FOR INSTITUTIONAL, HOSPITAL AND FARM PREMISE USE.

• **Virucidal**:** Avian Influenza A (H5N2), Influenza A₂ (Japan 305/57 Asian strain), Cytomegalovirus, Coronavirus, Polio Type 1 virus, Vaccinia, Canine parvovirus, Herpes simplex types 1/2 and 2/G (oral, ocular and genital), HIV-1 (AIDS Virus). See Sporidicin Bulletin NO. 303.

• **Bactericidal:** Staphylococcus aureus, Salmonella enterica, Pseudomonas aeruginosa, Methicillin resistant Staphylococcus aureus (MRSA), Vancomycin resistant Enterococcus faecium (VRE), E. coli, Streptococcus pyogenes, Streptococcus salivarius, Mycobacterium tuberculosis.

• **Fungicidal:** Chaetomium globosum, Aspergillus niger, Penicillium variable, Trichophyton mentagrophytes (Athlete's Foot fungus).

ACTIVE DEODORIZER*: Provides continuous residual bacteriostatic activity against odor-causing bacteria for up to 6 months, if not removed from surfaces including plastics, latex, vinyl, glass, treated wood, metal, glazed porcelain, glazed tile and paint.


SPORICIDIN DISINFECTANT TOWELETTES are germicidal when used as directed for disinfection but differ in composition from other Sporidicin products.

DIRECTIONS FOR USE: It is a violation of U.S. Federal Law to use this product in a manner inconsistent with its labeling.

DO NOT DILUTE OR MIX WITH OTHER CHEMICALS

PRECLEAN surfaces using Sporidicin Disinfectant Solution to remove soil and filth. Wipe dry with paper towel, cloth or sponge.

TO DISINFECT AND DEODORIZE items which cannot be immersed, such as electrical panels, cash registers, door knobs, handles, seats, dashboards and other "high touch" surfaces in and around farm equipment, animal and poultry housing, barns, kennels, breeding pens, hatcheries, trucks and vehicles: Thoroughly wet pre-cleaned surface with

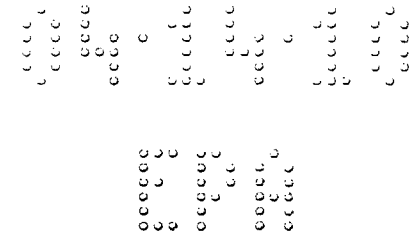

Sporidicin Disinfectant Towelettes and allow surface to remain wet for 10 minutes. Use as many towelettes as necessary for the treated area to remain wet for ten minutes at room temperature (20°C/68°F or above) to kill listed organisms.

TO KILL, CLEAN AND PREVENT MOLD AND MILDEW: Use procedure stated above. Repeat application if new growth appears.

TO CLEAN AIR DUCTS: Apply to hard, non-porous surfaces.

TO SPOT CLEAN, DEODORIZE AND REMOVE DEBRIS FROM CARPETS AND FABRICS: Blot or scrub area to be treated with Sporidicin Disinfectant Towelette to remove soiling. Wipe dry with clean cloth or towel.

FOR FOOD PROCESSING PLANTS: Follow pre-cleaning and disinfecting instructions outlined above to disinfect hard, non-porous surfaces and materials in food handling areas, including poultry and meat packing facilities, slaughter houses, etc. All surfaces and materials that come in contact with food must be washed with soap or detergent and rinsed with potable water prior to use.

PRECAUTIONARY STATEMENTS: HAZARD TO HUMANS AND DOMESTIC ANIMALS.

CAUTION: Causes moderate eye irritation. Avoid contact with eyes or clothing. Wash thoroughly with soap and water after handling.

FIRST AID: If in eyes, remove contact lenses. Hold eye open and slowly rinse with water. If irritation persists for more than 5 minutes continue rinsing and consult a physician or poison control center.

STORAGE AND DISPOSAL: Store at room temperature. Nonrefillable container. Do not reuse or refill this container. Discard empty container in trash or offer for recycling if available.

FOR HOSPITAL DISINFECTION: This product is not to be used as a terminal sterilant/high level disinfectant on any surface or instrument that (1) is introduced directly into the human body, either into or in contact with the bloodstream or normally sterile areas of the body, or (2) contacts intact mucous membranes but which does not ordinarily penetrate the blood barrier or otherwise enter normally sterile areas of the body. This product may be used to pre-clean or decontaminate critical or semi-critical medical devices prior to sterilization or high level disinfection.

Manufactured for
Contec, Inc. / Sporidicin Division
Spartanburg, SC 29303
1-800-424-3733 www.sporidicin.com

AF0020-03/10

NOTIFICATION
Date Reviewed: 4/16/10
Reviewed By: S. Grady

474