

34810-25

05/13/2010

1/6

UNITED STATES ENVIRONMENTAL PROTECTION AGENCY



UNITED STATES ENVIRONMENTAL PROTECTION AGENCY  
WASHINGTON, DC 20460

MAY 13 2010

OFFICE OF  
PREVENTION, PESTICIDE  
AND TOXIC SUBSTANCE

Mary Anne Auer  
**Wexford Labs, Inc.**  
325 Leffingwell Avenue  
Kirkwood, MO 63122

Subject: **Ready To Use Thymo-Cide**  
EPA Registration No. 34810-25  
Application Date: February 10, 2010  
Receipt Date: February 12, 2010

Dear Ms. Auer:

The following amendment, submitted in connection with registration under the Federal Insecticide, Fungicide, and Rodenticide Act (FIFRA), as amended, is accepted with conditions.

**Proposed Amendment:**

- Add Pandemic Influenza A H1N1
- List Deodorizing Claim separately
- Addition of optional marketing language that includes additional types of surfaces and use sites
- Option Signal Word
- Update label as per PR Notice 2007-4

**Conditions:**

Revise the label as follows:

- 1.) Delete the phrase "disease causing" on page 2 of your label
- 2.) Revise the Personal Protection Statement under HIV instructions on page 4 by deleting the word "or" and stating "and", also replace the word "should" with "must" to be in compliance with PRN 2005-1.
- 3.) Revise your Storage and Disposal statements to be in compliance with 2007-4 for the institutional/commercial use sites cited on the label by stating the following:

**Storage and Disposal**

Do not contaminate food, feed, or water by storage or disposal.

**Pesticide Storage:** Store at room temperature. Do not leave containers of this product open during storage. Store in original container in areas inaccessible to children and persons unfamiliar with its use.

**Pesticide Disposal:** Wastes resulting from the use of this product must be disposed of on site or at an approved waste disposal facility.

**Container Handling:** Proposed statements are acceptable

- 4.) Delete the phrase, containing botanical active ingredient found in nature, in the second paragraph under "Cleaning and Disinfection Instructions" on page 2


SYMBOL							
SURNAME							
DATE							

**General Comment:**

A stamped label is enclosed. Submit 1 copy of your final printed label before distributing or selling the product bearing the revised labeling.

Should you have any questions or comments concerning this letter, you may contact me by telephone at (703) 308-6416 or by e-mail at [campbell-mcfarlane.jacqueline@epa.gov](mailto:campbell-mcfarlane.jacqueline@epa.gov) or Stacey Grigsby by telephone at (703) 305-6440 or by e-mail at [grigsby.stacey@epa.gov](mailto:grigsby.stacey@epa.gov) during the hours of 8:00am to 4:00pm EST. When submitting information or data in response to this letter, a copy of this letter should accompany the submission to facilitate processing.

Sincerely,

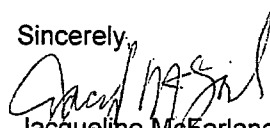

Jacqueline McFarlane  
 Product Manager (34)  
 Regulatory Management Branch II  
 Antimicrobials Division (7510P)

**Enclosure: EPA stamped accepted labeling**

CONCURRENCES							
SYMBOL							
SURNAME							
DATE							

{All text in brackets [xxx] is optional & may or may not be included on a final label.}  
{All text in braces {xxx} is administrative & will not appear on a final label.}

**READY TO USE THYMO-CIDE**

[ONE-STEP] GERMICIDAL DETERGENT  
[SERUM EFFECTIVE ONE-STEP GERMICIDAL DETERGENT]  
CLEANS • DISINFECTS • DEODORIZES

GERMICIDE [GERMICIDAL]  
[PSEUDOMONACIDAL]  
[BACTERICIDAL] [BACTERICIDE]  
TUBERCULOCIDE [TUBERCULOCIDAL]  
FUNGICIDE [FUNGICIDAL]  
VIRUCIDE\* [VIRUCIDAL\*] [with 5% Serum]

[Kills Pandemic 2009 H1N1 Influenza A Virus (formerly called swine flu)]  
[Hospital Disinfectant]  
[Instrument PreSoak]  
[Disinfectant Cleaner]  
[Decontaminant]  
[Herbal or Citrus Scent]

ACCEPTED  
with COMMENTS  
in EPA Letter Dated:

**MAY 13 2010**

Under the Federal Insecticide,  
Fungicide, and Rodenticide Act as  
amended, for the pesticide,  
registered under EPA Reg. No. **34810-25**

**ACTIVE INGREDIENT:**

Thymol ..... 0.051%

**INERT INGREDIENTS:** ..... 99.949%

[Includes detergents, perfume and water] 100.000%

**[CAUTION]**

**KEEP OUT OF REACH OF CHILDREN**

**[FIRST AID**

**FOR SKIN:** Take off contaminated clothing. Rinse skin immediately with plenty of water for 15-20 minutes. Call a Poison Control Center or doctor for treatment advice. Have the product container or label with you when calling a poison control center or doctor or going for treatment.]

**[PRECAUTIONARY STATEMENTS: HAZARD TO HUMANS AND DOMESTIC ANIMALS. CAUTION:** Harmful if absorbed through skin. Avoid contact with eyes, skin or clothing. Wash hands before eating, drinking, chewing gum, using tobacco or using the toilet.]

Wexford Labs, Inc.  
325 Leffingwell Avenue • Kirkwood, Missouri 63122

EPA REG. NO. 34810-25 • EPA EST. NO. 34810-MO-1  
Net Contents: \_\_\_\_ [oz] [gal] [ml]

**GERMICIDAL ACTIVITY:** (AOAC Use Dilution Confirmation in the presence of 5% serum at 20° C for 10 minutes.) Pseudomonas aeruginosa [Ps. Aeruginosa] (ATCC 15442), Staphylococcus aureus [S. aureus] (ATCC 6538), Salmonella enterica [S. enterica] (ATCC 10708)

**BROAD SPECTRUM DATA:** (AOAC Use Dilution Confirmation in the presence of 5% serum at 20° C for 10 minutes.) Methicillin resistant Staphylococcus aureus (MRSA)\*\* , Methicillin resistant Staphylococcus epidermidis (MRSE)\*\* Tetracycline, Novobiocin Streptomycin resistant Staphylococcus aureus (TNSR) (ATCC 27659), Gentamicin, Methacillin resistant Staphylococcus aureus (GRM) (ATCC 33592), Penicillin resistant Staphylococcus (PR) (ATCC 11631), Enterococcus Faecalis (VRE) (ATCC 51299).  
\*\* Clinical Isolates

**TUBERCULOCIDAL ACTIVITY:** Tuberculocidal in the presence of 5% serum at 20° C for 10 minutes. AOAC Mycobacterium bovis [M. tuberculosis] [TB].

**FUNGICIDAL ACTIVITY:** AOAC Trichophyton Mentagrophytes [T. Mentagrophytes] (ATCC 9533) in the presence of 5% serum at 20° C for 10 minutes.

**VIRUCIDAL ACTIVITY:** \*Influenza A3 (Hong Kong), Herpes Simplex Type 2 (VR 734), Vaccinia Virus, Human Rotavirus, and Respiratory Syncytial Virus (VR 26) in the presence of 5% serum at 20° C for 10 minutes. Human immunodeficiency virus Type 1 (HIV-1) (AIDS) in the presence of 10% serum on hard, non-porous surfaces at 20° C for 10 minutes. This product has demonstrated effectiveness against influenza A3 virus and is expected to inactivate all influenza A viruses including Pandemic 2009 H1N1 (formerly called swine flu).

**DIRECTIONS FOR USE**

It is a violation of Federal Law to use the product in a manner inconsistent with its labeling.

**CLEANING AND DISINFECTING DIRECTIONS:** Gross filth and heavy soil must be removed before applying cleaning solution. Apply to surfaces by cloth, sponge, brush, coarse spray or by immersing equipment using normal cleaning methods. All hard surfaces must be wet thoroughly, and remain wet for at least 10 minutes and allowed to air dry. Thoroughly rinse all wetted and cleaned food contact surfaces with potable water. As the solution becomes dirty discard and replace with fresh solution. This is a complete product. DO NOT MIX WITH OTHER CHEMICALS.

Ready To Use Thymo-Cide is a [synthetic] germicidal detergent [containing a botanical active ingredient found in nature] [with broad spectrum kill of gram positive and gram negative microorganisms] formulated to clean, disinfect and deodorize in a simple one-step process. [It is an excellent [one-step] cleaner [disinfectant] for working surfaces, including stainless steel, chrome, glass and vinyl]. [It effectively controls odors produced by odor causing bacteria and leaves a fresh citrus scent.]

Use Ready To Use Thymo-Cide for all disinfection of hard, non-porous surfaces (floors, walls, tables, etc.) in Healthcare Facilities including [Operating Rooms], [Intensive Care],

[Nurseries], [Emergency Areas], [Dental Operatories], [Dental and Medical Offices], [Oral Healthcare Facilities], [Oral Surgery Centers], [Dental Schools], [Schools], [Hospitals], [Nursing Homes], [Prisons], [Dental Labs], [by Police and EMS], as well as other areas/facilities concerned with the hazards of cross contamination.

[Use for one-step cleaning, disinfecting and deodorizing of devices and surfaces potentially contaminated with disease causing microorganisms and other items made of stainless steel, chrome, glass, non-porous vinyl, glazed porcelain, metal, plastics, Plexiglas®, glazed tile, washable painted or varnished surfaces, linoleum, and rubber.]

[When used as directed, Ready To Use Thymo-Cide effectively cleans and disinfects hard, non-porous environmental surfaces such as: ambulance patient care surfaces, anesthesia equipment surfaces, baby/infant care surfaces, backboards, bathing units, bathrooms, bed railings, blood pressure devices, cabinets, chairs, changing tables, computer surfaces, countertops, CT equipment/surfaces, curing lights, dental equipment surfaces/devices, doorknobs, floors, garbage cans, highchairs, infant care equipment, keyboards, laboratory equipment/surfaces, lights, light lens covers, lumens, medical equipment surfaces/devices, microscopes, MRI equipment surfaces, neonatal/nursery surfaces, obstetrical equipment surfaces, operating room lights/tables, ophthalmic equipment surfaces/devices, optical wear (excluding contact lenses), oxygen hoods, patient bed surfaces, patient stretchers, patient transport surfaces, physical therapy equipment, play-care equipment, play mats, radiology equipment surfaces, refrigerators (exterior surface), respirators, respiratory care equipment surfaces, sinks, soaking containers, stalls, shower floors/walls, skin care equipment/surfaces, surgical equipment surfaces, surgical microscopes, tables, telephones, toilets (rims, seats and exterior surfaces), toys, transportation equipment surfaces, trays, walls, work stations, x-ray equipment/surfaces, and other areas/facilities concerned with the hazards of cross contamination.]

**ODOR CONTROL:** It effectively controls odors by killing odor-causing bacteria while simultaneously chemically neutralizing their odors and leaving a pleasant fragrance. Follow "CLEANING AND DISINFECTING DIRECTIONS".

**DECONTAMINATION OF CRITICAL OR SEMI-CRITICAL DEVICES PRIOR TO TERMINAL STERILIZATION/HIGH LEVEL DISINFECTION:** This product is not to be used as a terminal sterilant/high level disinfectant on any surface or instrument that (1) is introduced directly into the human body, either into or in contact with the bloodstream or normally sterile areas of the body, or (2) contacts intact mucous membranes but which does not ordinarily penetrate the blood barrier or otherwise enter normally sterile areas of the body. This product may be used to preclean or decontaminate critical or semi-critical medical devices prior to sterilization or high-level disinfection.

**\*READY TO USE THYMO-CIDE KILLS HIV ON PRE-CLEANED ENVIRONMENTAL SURFACES/OBJECTS PREVIOUSLY SOILED WITH BLOOD/BODY FLUIDS** in all settings in which there is an expected likelihood of soiling of hard, non-porous surfaces/objects with blood or body fluids, and in which the surfaces/objects likely to be soiled with blood or body fluids can be associated with the potential for transmission of human immunodeficiency virus type 1 (HIV-1) (associated with AIDS).

**SPECIAL INSTRUCTIONS FOR CLEANING AND DECONTAMINATING SURFACES/OBJECTS SOILED WITH BLOOD AND BODY FLUIDS WHICH MAY CONTAIN HIV-1 (AIDS virus).** **PERSONAL PROTECTION:** Wear disposable rubber or latex gloves, gowns, face masks or eye coverings as appropriate, when performing manual cleaning of these surfaces. **CONTACT TIME:** Leave surfaces wet for 1 minute. This contact time will not control other common types of viruses and bacteria. **CLEANING PROCEDURE:** The solution may be used to remove and clean blood and other body fluids from surfaces and objects before disinfection. For lightly soiled areas wipe off visual soil/clean with solution, and disinfect surface with a second application of the solution. For large spills the contaminated area should be flooded with the solution before cleaning, then decontaminated with fresh solution. **DISPOSAL OF INFECTIOUS MATERIAL:** Blood and other body fluids should be autoclaved and disposed of according to Federal State and local regulations for infectious waste disposal.

**STORAGE AND DISPOSAL:**

Store at room temperature. Do not leave containers of this product open during storage. Store in original container in areas inaccessible to children.

Unused product may be discarded by pouring down the drain and flushing drain with water.

{For non-refillable container:} Nonrefillable container. Do not reuse or refill this container. Clean container promptly after emptying. Triple rinse as follows: Fill container 1/4 full with water and reclose the container. Agitate vigorously and dispose of rinsate. Repeat two more times. [Wrap container and put in trash or] offer for recycling, if available.

{For refillable container:} Refillable container. Refill this container with pesticide only. Do not reuse this container for any other purpose. Cleaning the container before final disposal is the responsibility of the person disposing of the container. Cleaning before refilling is the responsibility of the refiller. To clean the container before final disposal, triple rinse container (or equivalent) promptly after emptying. Offer for recycling, if available.