

US ENVIRONMENTAL PROTECTION AGENCY OFFICE OF PESTICIDES PROGRAMS REGISTRATION DIVISION (TS-767) WASHINGTON, DC 20460 NOTICE OF PESTICIDE: <input checked="" type="checkbox"/> REGISTRATION <input type="checkbox"/> REREGISTRATION <i>(Under the Federal Insecticide, Fungicide, and Rodenticide Act, as amended)</i>	EPA REGISTRATION NO. 39860-2	DATE OF ISSUANCE DEC 17 1985
	TERM OF ISSUANCE	
	NAME OF PESTICIDE PRODUCT Chemcide	
NAME AND ADDRESS OF REGISTRANT (Include ZIP code) <div style="border: 1px solid black; padding: 10px; width: fit-content; margin: auto;"> <p>Chemrion Laboratories, Inc. 5211 Telegraph Road Los Angeles, CA 90022</p> </div>		
NOTE: Changes in labeling formula differing in substance from that accepted in connection with this registration must be submitted to and accepted by the Registration Division prior to use of the label in commerce. In any correspondence on this product always refer to the above U.S. EPA registration number.		
On the basis of information furnished by the registrant, the above named pesticide is hereby Registered/Reregistered under the Federal Insecticide, Fungicide, and Rodenticide Act. A copy of the labeling accepted in connection with this Registration/Reregistration is returned herewith. Registration is in no way to be construed as an indorsement or approval of this product by this Agency. In order to protect health and the environment, the Administrator, on his motion, may at any time suspend or cancel the registration of a pesticide in accordance with the Act. The acceptance of any name in connection with the registration of a product under this Act is not to be construed as giving the registrant a right to exclusive use of the name or to its use if it has been covered by others. This product is conditionally registered in accordance with FIFRA sec. 3(c)(7)(A) provided that you: <ol style="list-style-type: none"> 1. Submit and/or cite all data required for registration/reregistration of your product under FIFRA sec. 3(c)(5) when the Agency requires all registrants of similar products to submit such data. 2. Make the labeling changes listed below before you release the product for shipment: <ol style="list-style-type: none"> a. Add the phrase "EPA Registration No. 39860-2." b. Revise the statement: <p style="margin-left: 40px;">Hazardous to Humans . . .</p> <p style="margin-left: 40px;">to read:</p> <p style="margin-left: 40px;">Hazards to Humans</p> 		
<input checked="" type="checkbox"/> ATTACHMENT IS APPLICABLE		
SIGNATURE OF APPROVING OFFICIAL		DATE

c. Revise the heading:

Precaution

to read:

Precautionary Statements

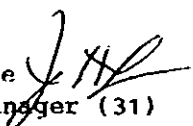
d. Delete the statement:

Brand of Benzalkonium Chloride

3. Submit five (5) copies of your final printed labeling before you release the product for shipment. Refer to the A-79 Enclosure for a further description of final printed labeling.

If these conditions are not complied with, the registration will be subject to cancellation in accordance with FIFRA sec. 6(e). Your release for shipment of the product constitutes acceptance of these conditions.

A stamped copy of the label is enclosed for your records.

John H. Lee 
Product Manager (31)
Disinfectants Branch
Registration Division (TS-767C)

Enclosure

FRONT PANEL

Benzalkonium Chloride
Aqueous Solution 1:750
For Disinfection of Dental and
Surgical Instruments

Active Ingredient:
Alkyl (40% C₁₂, 50% C₁₄, 10% C₁₆) Dimethyl
Benzyl Ammonium Chloride 0.133%
Inert Ingredients 99.867%
CONTAINS SODIUM NITRITE

KEEP OUT OF REACH OF CHILDREN
CAUTION

See side panel for additional precautionary statement.

E. P. A. REG. NO.

Mfg. by
CHEM-RICH LABORATORIES, INC.
Los Angeles, CA 90022

Net Contents: One Gallon

E. P. A. EST. NO. 39860-CA-001

ACCEPTED
FOR COMMENTS

DEC 31 1985

Under the Federal Insecticide,
Fungicide, and Rodenticide Act
as amended, for the pesticide
registered under EPA Reg. No.

39860-2

SIDE PANEL

An efficient Chemical Germicide for the Disinfection of
Dental and Surgical Instruments.

Non Rusting Stainless Stable

NOTE

This will not rust or corrode surgical steel, or other
metals. May be used for disinfection of metal, plastic,
rubber or glass instruments and apparatus including
catheters, gloves, hypodermic equipment and tubing.

Not to be relied to destroy bacterial spores, Mycobacterium
tuberculosis or the etiological agent of viral hepatitis.

PRECAUTION
Hazardous to Human and Domestic Animals

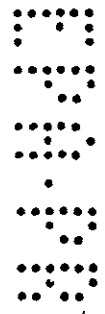
DIRECTION FOR USE

It is a violation of Federal Law to use this product
inconsistent with its labeling.
The solution should be used full strength. Clean instruments
well and rinse. Submerge article in solution for at least
10 minutes. Flush and fill cleansed lumen of any hollow
instrument. Instrument may be used immediately without further
rinsing or drying. When rinsing is necessary, use sterile
water and aseptic procedure for critical items. Use a fresh
solution daily or when solution becomes diluted or soiled.

For destruction of spore forming bacteria, steam autoclaving
is mandatory, after which instruments or needles may be stored
in a sterile covered pan in this solution until used.

DISPOSAL

Do not reuse empty container. Rinse empty container thoroughly
with water and discard in trash.



c. Revise the heading:

Precaution

to read:

Precautionary Statements

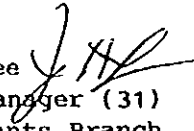
d. Delete the statement:

Brand of Benzalkonium Chloride

3. Submit five (5) copies of your final printed labeling before you release the product for shipment. Refer to the A-79 Enclosure for a further description of final printed labeling.

If these conditions are not complied with, the registration will be subject to cancellation in accordance with FIFRA sec. 6(e). Your release for shipment of the product constitutes acceptance of these conditions.

A stamped copy of the label is enclosed for your records.

John H. Lee 
Product Manager (31)
Disinfectants Branch
Registration Division (TS-767C)

Enclosure

FRONT PANEL

Benzalkonium Chloride
Aqueous Solution 1:750
For Disinfection of Dental and
Surgical Instruments

Active Ingredient:
Alkyl (40% C₁₂, 50% C₁₄, 10% C₁₆) Dimethyl
Benzyl Ammonium Chloride 0.133%
Inert Ingredients 99.867%
CONTAINS SODIUM NITRITE

KEEP OUT OF REACH OF CHILDREN
CAUTION

See side panel for additional precautionary statement.

E. P. A. REG. NO.
Mfg. by
CHEMRICH LABORATORIES, INC.
Los Angeles, CA 90022
Net Contents: One Gallon

ACCEPTED
WITH COMMENTS
DEC 3 1985
Under the Federal Insecticide,
Fungicide, and Rodenticide Act
as amended, for the pesticide
registered under EPA Reg. No.
39860-2

E. P. A. EST. NO. 39860-CA-001

SIDE PANEL

An efficient Chemical Germicide for the Disinfection of
Dental and Surgical Instruments.

Non Rusting Stainless Stable

NOTE

This will not rust or corrode surgical steel, or other
metals. May be used for disinfection of metal, plastic,
rubber or glass instruments and apparatus including
catheters, gloves, hypodermic equipment and tubing.

Not to be relied to destroy bacterial spores, Mycobacterium
tuberculosis or the etiological agent of viral hepatitis.

PRECAUTION
Hazardous to Human and Domestic Animals

DIRECTION FOR USE

It is a violation of Federal Law to use this product
inconsistent with its labeling.
The solution should be used fullstrength. Clean instruments
well and rinse. Submerge article in solution for at least
10 minutes. Flush and fill cleansed lumen of any hollow
instrument. Instrument may be used immediately without further
rinsing or drying. When rinsing is necessary, use sterile
water and aseptic procedure for critical items. Use a fresh
solution daily or when solution becomes diluted or soiled.

For destruction of spore forming bacteria, steam autoclaving
is mandatory, after which instruments or needles may be stored
in a sterile covered pan in this solution until used.

DISPOSAL

Do not reuse empty container. Rinse empty container thoroughly
with water and discard in trash.

