

**DIRECTIONS FOR USE**

It is a violation of Federal Law to use this product in a manner inconsistent with its labeling.

Use at 1:100 dilution (approximately 1 1/4 ounces per gallon of water). Sanco-Phene will provide 850 ppm phenols to this suggested use solution. Sanco-Phene should be applied by wiping, mopping, or scrubbing surfaces to wet thoroughly with the recommended use dilution. Allow surfaces to remain wet for 10 minutes to ensure effectiveness. Air dry. For areas of gross filth or heavy soil, a preliminary cleaning should be done using 1 1/4 ounces of Sanco-Phene per gallon of water. Then proceed as above.

This product is not to be used as a terminal sterilant/high level disinfectant on any surface or instrument that (1) is introduced directly into the human body, either into or in contact with the bloodstream or normally sterile areas of the body, or (2) contacts intact mucous membranes but which does not ordinarily penetrate the blood barrier or otherwise enter normally sterile areas of the body. This product may be used to pre-clean or decontaminate critical or semi-critical medical devices prior to sterilization or high level disinfection.

**STORAGE AND DISPOSAL**

**Storage:** Store in original container, tightly closed and in a safe place. Do not use or store near heat or open flame.

**Pesticide Disposal:** Pesticide wastes are acutely hazardous. Improper disposal of excess pesticide, spray-mixture or rinsate is a violation of Federal Law. If these wastes cannot be disposed of by use according to label instructions, contact your State Pesticide or Environmental Control Agency or the Hazardous Waste representative at the nearest EPA Regional Office for guidance. Do not contaminate water, food or feed by storage or disposal.

**Container Disposal:** Triple-rinse (or equivalent) then offer for recycling or reconditioning or puncture and dispose of in a sanitary landfill, or incineration, or, if allowed by state and local authorities, by burning. If burned, stay out of smoke.

**SANCO-PHENE  
CLEANER DISINFECTANT DEODORANT**

**ACTIVE INGREDIENTS:**

- o-Benzyl-p-chlorophenol ..... 4.00%
- o-Phenylphenol ..... 3.00%
- p-Tertiary amylphenol ..... 1.50%

**INERT INGREDIENTS\*** ..... 91.50%  
100.00%

(\* Contains cleaning and water softening agents)

**KEEP OUT OF REACH OF CHILDREN**

**DANGER**

See side panel for additional precautionary statements.

**STATEMENT OF PRACTICAL TREATMENT**

If on skin - wash with plenty of soap and water. Get medical attention if irritation persists.

If in eyes - immediately flush eyes with plenty of water for at least 15 minutes. Get medical attention.

If swallowed - promptly drink large quantities of water. Avoid alcohol. Call a physician immediately.

**NOTE TO PHYSICIAN**

Probable mucosal damage may contraindicate the use of gastric lavage. Measures against circulatory shock, respiratory depression and convulsion may be needed.

EPA REG. NO. 49403-5  
EPA EST. NO. 49403-EN-1

**NET CONTENTS:**

SANCO-PHENE is a concentrated phosphate-free, germicidal cleaner intended for cleaning and disinfection. It is designed for cleaning floors, walls, bed frames, furniture, and bathroom fixtures. It is free-rinsing and does not leave a film. Consequently, it is recommended for the cleaning of conductive flooring in operating rooms.

SANCO-PHENE disinfectant has been tested according to the AOAC Use Dilution Method, Official Methods of Analysis and proven effective against:

- Staphylococcus aureus
- Salmonella choleraesuis
- Pseudomonas aeruginosa

**PRECAUTIONARY STATEMENTS**

**HAZARDS TO HUMANS AND DOMESTIC ANIMALS**


**DANGER.** Corrosive. Causes eye damage and severe skin irritation. Do not get in eyes, on skin or on clothing. Protect eyes and skin when handling concentrate by wearing goggles and rubber gloves. Harmful or fatal if swallowed. Avoid contamination of food. Wash thoroughly with soap and water after handling and before eating or smoking. Remove contaminated clothing and wash before reuse.

**PHYSICAL AND CHEMICAL HAZARDS**

Do not use or store near heat or open flame.

**MANUFACTURED BY:**

Nipa Laboratories, Inc.  
3411 Silverside Road  
104 Hagley Building  
Wilmington, DE 19810

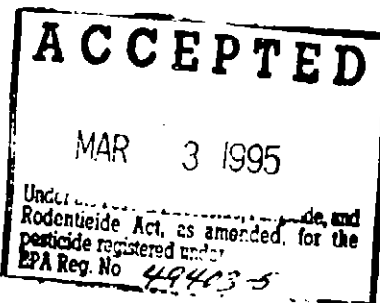

**BEST AVAILABLE COPY**