

UNITED STATES ENVIRONMENTAL PROTECTION AGENCY

Andersen Sterilizers, Inc.
3154 Caroline Dr.
Haw River, NC 27258

MAY 10 2000

Subject: Anprolene AN-79
EPA Registration No. 69340-2

EOGas AN-1004/5/6
EPA Registration No. 69340-4

Amendment Date: 2/3/00
EPA Receive Date: 2/10/00

Attention: Mr. Juergen Ballo

The amendments, referenced above, submitted in connection with registration under FIFRA, as amended, to add new packaging configurations for EOGas Commercial Refill Packs are acceptable provided you address the following recommendations before you release this product for shipment.

Label Recommendations:

- a) Under the "Precautionary Statements," you will need to add the statements, "Do not breathe vapor. Do not swallow," after the "Harmful if inhaled" statement.
- b) The EPA Registration Number for this product must be on the label.


Manual Recommendation

The heading, "Instructions for Use," must be revised to state "Directions for Use."

Stamped labels are enclosed for your records. You must submit a final printed copy of the label and manual prior to the release of the products for shipment. If you have any questions regarding this letter, please contact Jacqueline Campbell at (703) 308-6416.

CONCURRENCES							
SYMBOL	7570C						
SURNAME	Campbell						
DATE	5/10/00						

Sincerely,

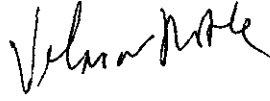A handwritten signature in cursive script that reads "Velma Noble".

Velma Noble
Product Manager (31)
Regulatory Management Branch I
Antimicrobials Division (7510C)

ACCEPTED
with COMMENTS
in EPA Letter Dated:

MAY 10 2000

Y. Registrado como Fungicida,
Herbicida y Bactericida Act as
Fungicide, Herbicide, Bactericide,
registered under EPA Reg. No. 69340-4




AN 2003

KEEP OUT OF REACH OF CHILDREN DANGER

IMPORTANT: For Professional Use Only
Users must follow the requirements of the OSHA occupational exposure standard for Ethylene oxide (29 CFR 1910.1047)

Active ingredient: Ethylene oxide	89.0%
Inert ingredients:	11.0%
TOTAL	100.0%

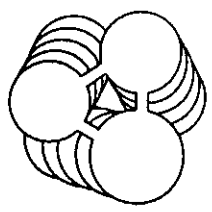

EO Gas™

INDUSTRIAL REFILL PACK

ONE PACK

CONTENTS: 50 Cartridges. Each AN2003 cartridge contains 0.09 av. oz. (3.0 g) EO Gas sterilant.

ANDERSEN STERILIZERS, INC.

Manufactured and Distributed by:
ANDERSEN STERILIZERS, INC.
3154 Caroline Drive • Haw River, NC 27258 USA
336-376-8622

Distributed by:
ANDERSEN PRODUCTS, INC.
3202 Caroline Drive • Haw River, NC 27258 USA
1-800-523-1276 • 336-376-3000
www.anpro.com • mailbox@anpro.com

European Distributor:
H. W. ANDERSEN PRODUCTS, LTD.
Clacton-On-Sea • Essex CO15 4XA UK

EPA Registration No. XXXXX-X
EPA Establishment No. 9417-NC-001



PRECAUTIONARY STATEMENTS

HAZARDS TO HUMANS AND DOMESTIC ANIMALS.
DANGER! CANCER HAZARD AND REPRODUCTIVE HAZARD.
CAUSES IRREVERSIBLE EYE DAMAGE AND SKIN BURNS.
HARMFUL IF INHALED. MAY CAUSE NERVOUS SYSTEM DAMAGE. DO NOT GET IN EYES, SKIN, OR ON CLOTHING.
STORE AND USE WITH ADEQUATE VENTILATION.

See user's manual for additional precautionary statements.

STATEMENT OF PRACTICAL TREATMENT

In all cases of exposure, get medical attention immediately. Take person to a doctor or emergency treatment facility.

IF IN EYES: Hold eyelids open and flush with a steady, gentle stream of water for at least 15 minutes.

IF INHALED: Get exposed person to fresh air. Keep warm. Make sure person can breathe freely. If breathing is difficult, give oxygen. Call a physician even if no symptoms are present. Keep under medical observation. Symptoms may be delayed. If breathing has stopped, give artificial respiration. Do not give anything by mouth to an unconscious person.

IF ON SKIN: Immediately wash skin with plenty of water while removing contaminated clothing and shoes. Call a physician. Aerate, wash, or clean contaminated clothing and discard leather goods.

IF SWALLOWED: Drink at least two glasses of water. Do not induce vomiting. Do not give anything by mouth to an unconscious victim. Call a physician.

NOTE TO PHYSICIAN

Ethylene oxide is a gas. Skin exposure from contact with fabric, rubber or plastic containing residual ethylene oxide will commonly result in skin irritation with extensive blister formation. At high concentrations severe conjunctivitis can occur. Irritation of the respiratory tract may occur, but without acute lung edema. Symptoms of systemic intoxication are headache, nausea, vomiting, incoordination, and cardiac irregularities. Treatment is symptomatic.

ENVIRONMENTAL HAZARDS

Do not discharge effluent containing this product into lakes, streams, ponds, estuaries, oceans, or other waters unless in accordance with the requirements of a National Discharge Eliminations System (NDSES) permit and the permitting authority has been notified in writing prior to discharge. Do not discharge effluent containing this product to sewer systems without previously notifying the local sewage treatment plant authority. For guidance contact your State Water Board or Regional Office of the EPA.

PHYSICAL OR CHEMICAL HAZARDS

Ethylene oxide gas is extremely flammable. Do not use near flame, electrical sparks, or hot surfaces. Ground all equipment to prevent static sparks.

DIRECTIONS FOR USE

It is a violation of Federal law to use this product in a manner inconsistent with its labeling. It is the responsibility of the employer of any person engaged in the handling or application of this product to follow the requirements of 29 CFR 1910.1047. **To be used only in an EO Gas sterilizer for hospital, medical, and veterinary sterilization.** Read package insert for complete sterilization instructions and additional precautions.

STORAGE AND DISPOSAL

Do not contaminate water, food, or feed by storage or disposal.
STORAGE: Store in a cool place and out of direct sunlight.

PESTICIDE DISPOSAL: Unwanted or expired EO Gas cartridges should be returned to the manufacturer for disposal. Contact Andersen Products, Inc. (tel. no. 336-376-3000) for instructions. If unwanted or expired EO Gas cartridges cannot be disposed of according to label instructions, contact your State Pesticide or Environmental Control Agency, or the Hazardous Waste Representative at the nearest EPA Regional Office for guidance.

CONTAINER DISPOSAL: Do not reuse empty box, used cartridge, or used liner bags. Wrap box, cartridge, and liner bags in paper and discard in the trash.



EO Gas® is a trademark of Andersen Products, Inc.
EPA Reg. No. 69340-4 and 5,150,700.

P/N XXXX 012800

AN2003 Kit 012800

1

2/3/00, 1:29 PM

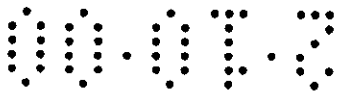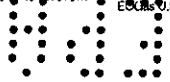

3/10



MADE IN
 U.S.A.
 in the USA Under Patent

MAY 10 2000
 Under the Federal Insecticide,
 Fungicide, and Rodenticide Act
 amended, for the pesticide,
 registered under EPA Reg. No. 69340-4

**AN
 2003**

EOGas™

REMOVE TAPE AND TRIGGER GUARD, PRESS BUTTON TO ACTIVATE CARTRIDGE

Use only in an Andersen EO Gas Sterilizer and according to manufacturer's instructions. See outer box and package insert for precautions and sterilization instructions.

Active Ingredient: Ethylene Oxide 89%
 Other Ingredient 11%
 Total 100%
 Net Contents: 25 Cartridges, 0.09 av. oz. (3.00 g)

LOT NO.

 USE
 BEFORE


DANGER
KEEP OUT OF REACH OF CHILDREN
CONTENTS FLAMMABLE
DO NOT OPEN THIS PACKAGE
DISPOSE OF SPENT EO Gas CARTRIDGE IN TRASH

CE 0413

Manufactured by: Andersen Sterilizers, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
 Distributed by: Andersen Products, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
 Distribué en Europe par: H. W. Andersen Products Ltd. • Clacton-On-Sea • Essex CO15 4XA U.K.
 EPA Registration No. XXXXX-X • EPA Est. No. 9417-NC-001 • U.S. Patent No. 4,937,046 & 5,160,700

PN XXXX 012800

REMOVE TAPE AND TRIGGER GUARD, PRESS BUTTON TO ACTIVATE CARTRIDGE

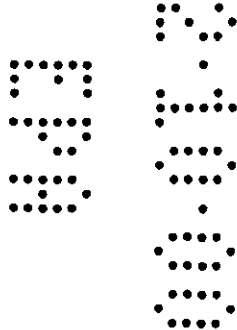

4/10

AN 2005

ONE PACK

CONTENTS: 50 Cartridges. Each AN2005 cartridge contains 0.18 av. oz. (5.0 g) EOGas sterilant.

ANDERSEN STERILIZERS, INC.

Manufactured and Distributed by:
ANDERSEN STERILIZERS, INC.
3154 Caroline Drive • Haw River, NC 27258 USA
336-376-8622

Distributed by:
ANDERSEN PRODUCTS, INC.
3202 Caroline Drive • Haw River, NC 27258 USA
1-800-523-1276 • 336-376-3000
www.anpro.com • mailbox@anpro.com

European Distributor:
H. W. ANDERSEN PRODUCTS, LTD.
Clacton-On-Sea • Essex CO15 4XA UK

EPA Registration No. XXXXX-X
EPA Establishment No. 9417-NC-001

KEEP OUT OF REACH OF CHILDREN

DANGER

IMPORTANT: For Professional Use Only

Users must follow the requirements of the OSHA occupational exposure standard for Ethylene oxide (29 CFR 1910.1047)

Active ingredient: Ethylene oxide 90.0%
Inert ingredients: 10.0%
TOTAL 100.0%

PRECAUTIONARY STATEMENTS

**HAZARDS TO HUMANS AND DOMESTIC ANIMALS.
DANGER! CANCER HAZARD AND REPRODUCTIVE HAZARD.
CAUSES IRREVERSIBLE EYE DAMAGE AND SKIN BURNS.
HARMFUL IF INHALED. MAY CAUSE NERVOUS SYSTEM
DAMAGE. DO NOT GET IN EYES, SKIN, OR ON CLOTHING.
STORE AND USE WITH ADEQUATE VENTILATION.**

See user's manual for additional precautionary statements.


STATEMENT OF PRACTICAL TREATMENT

In all cases of exposure, get medical attention immediately. Take person to a doctor or emergency treatment facility.

IF IN EYES: Hold eyelids open and flush with a steady, gentle stream of water for at least 15 minutes.

IF INHALED: Get exposed person to fresh air. Keep warm. Make sure person can breathe freely. If breathing is difficult, give oxygen. Call a physician even if no symptoms are present. Keep under medical observation. Symptoms may be delayed. If breathing has stopped, give artificial respiration. Do not give anything by mouth to an unconscious person.

IF ON SKIN: Immediately wash skin with plenty of water while removing contaminated clothing and shoes. Call a physician. Aerate, wash, or clean contaminated clothing and discard leather goods.

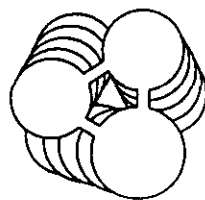

EOGas™

INDUSTRIAL REFILL PACK

IF SWALLOWED: Drink at least two glasses of water. Do not induce vomiting. Do not give anything by mouth to an unconscious victim. Call a physician.

NOTE TO PHYSICIAN

Ethylene oxide is a gas. Skin exposure from contact with fabric, rubber or plastic containing residual ethylene oxide will commonly result in skin irritation with extensive blister formation. At high concentrations severe conjunctivitis can occur. Irritation of the respiratory tract may occur, but without acute lung edema. Symptoms of systemic intoxication are headache, nausea, vomiting, incoordination, and cardiac irregularities. Treatment is symptomatic.

ENVIRONMENTAL HAZARDS

Do not discharge effluent containing this product into lakes, streams, ponds, estuaries, oceans, or other waters unless in accordance with the requirements of a National Discharge Eliminations System (NDPES) permit and the permitting authority has been notified in writing prior to discharge. Do not discharge effluent containing this product to sewer systems without previously notifying the local sewage treatment plant authority. For guidance contact your State Water Board or Regional Office of the EPA.

PHYSICAL OR CHEMICAL HAZARDS

Ethylene oxide gas is extremely flammable. Do not use near flame, electrical sparks, or hot surfaces. Ground all equipment to prevent static sparks.

DIRECTIONS FOR USE

It is a violation of Federal law to use this product in a manner inconsistent with its labeling. It is the responsibility of the employer of any person engaged in the handling or application of this product to follow the requirements of 29 CFR 1910.1047. To be used only in an EOGas sterilizer for hospital, medical, and veterinary sterilization. Read package insert for complete sterilization instructions and additional precautions.

STORAGE AND DISPOSAL

Do not contaminate water, food, or feed by storage or disposal.

STORAGE: Store in a cool place and out of direct sunlight.

PESTICIDE DISPOSAL: Unwanted or expired EOGas cartridges should be returned to the manufacturer for disposal. Contact Andersen Products, Inc. (tel. no. 336-376-3000) for instructions. If unwanted or expired EOGas cartridges cannot be disposed of according to label instructions, contact your State Pesticide or Environmental Control Agency, or the Hazardous Waste Representative at the nearest EPA Regional Office for guidance.

CONTAINER DISPOSAL: Do not reuse empty box, used cartridge, or used liner bags. Wrap box, cartridge, and liner bags in paper and discard in the trash.

EOGas® is a trademark of Andersen Products, Inc.
E.O. 12801 S. Pat. no. 4,937,046 and 5,160,700.

P/N XXXX 012800

AN2005 Kit 012800

CE 0413

ACCEPTED
with COMMERCIAL
in EPA Letter # 94-17-NC-001
MAY 10 2000

2/3/00, 11:27 AM

Under the Federal Insecticide,
Fungicide, and Rodenticide Act as
amended, for the pesticide,
registered under EPA Reg. No. 69340-4

5/10



**AN
2005**

EOGas™

Use only in an Andersen EOGas Sterilizer and according to manufacturer's instructions. See outer box and package insert for precautions and sterilization instructions.

Active Ingredient: Ethylene Oxide 90%
Other Ingredient 10%
Total.....100%
Net Contents: 0.18 av. oz. (5.00 g)

LOT NO.

USE
BEFORE

DANGER
KEEP OUT OF REACH OF CHILDREN
CONTENTS FLAMMABLE
DO NOT OPEN THIS PACKAGE
DISPOSE OF SPENT EOGas CARTRIDGE IN TRASH



Manufactured by: Andersen Sterilizers, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
Distributed by: Andersen Products, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
Distribué en Europe par: H. W. Andersen Products Ltd. • Clacton-On-Sea • Essex CO15 4XA U.K.
EPA Registration No. XXXX-X • EPA Est. No. 9417-NC-001 • U.S. Patent No. 4,937,046 & 5,160,700

PN XXXX 012800

REMOVE TAPE AND TRIGGER GUARD, PRESS BUTTON TO ACTIVATE CARTRIDGE

ACCEPTED
with 4 Dev. 4. 19. 20
in EPA Label 012800


MAY 10 2000

General Insecticides,
Insecticide Act 23
Pesticide,
EPA Reg. No. 169340-4

AN2005 Label 012800

1

2/3/00, 11:29 AM

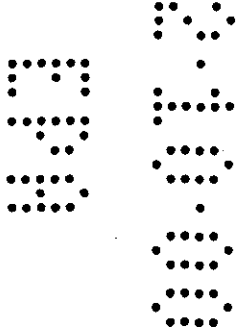

6/10

**AN
2011**


KEEP OUT OF REACH OF CHILDREN

DANGER

IMPORTANT: For Professional Use Only

Users must follow the requirements of the OSHA occupational exposure standard for Ethylene oxide (29 CFR 1910.1047)

Active Ingredient: Ethylene oxide 96.0%
Inert ingredients: 4.0%
TOTAL 100.0%

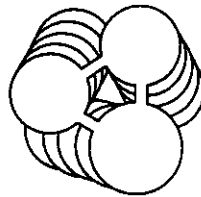

EOGas™
INDUSTRIAL REFILL PACK

ONE PACK

CONTENTS: 50 Cartridges. Each AN2011 cartridge contains 0.39 av. oz. (11.0 g) EOGas sterilant.

ANDERSEN STERILIZERS, INC.

Manufactured and Distributed by:

ANDERSEN STERILIZERS, INC.
3154 Caroline Drive • Haw River, NC 27258 USA
336-376-8622

Distributed by:

ANDERSEN PRODUCTS, INC.
3202 Caroline Drive • Haw River, NC 27258 USA
1-800-523-1276 • 336-376-3000
www.anpro.com • mailbox@anpro.com

European Distributor:

H. W. ANDERSEN PRODUCTS, LTD.
Clacton-On-Sea • Essex CO15 4XA UK

EPA Registration No. XXXXX-X
EPA Establishment No. 9417-NC-001



AN2011 Kit 012800

**ACCEPTED
with COMMENTS
in EPA Letter Dated:**

MAY 10 2000

Under the Federal Insecticide,
Fungicide, and Rodenticide Act as
amended, for the pesticide,
registered under EPA Reg. No. 69340-4

IF SWALLOWED: Drink at least two glasses of water. Do not induce vomiting. Do not give anything by mouth to an unconscious victim. Call a physician.

NOTE TO PHYSICIAN

Ethylene oxide is a gas. Skin exposure from contact with fabric, rubber or plastic containing residual ethylene oxide will commonly result in skin irritation with extensive blister formation. At high concentrations severe conjunctivitis can occur. Irritation of the respiratory tract may occur, but without acute lung edema. Symptoms of systemic intoxication are headache, nausea, vomiting, incoordination, and cardiac irregularities. Treatment is symptomatic.

ENVIRONMENTAL HAZARDS

Do not discharge effluent containing this product into lakes, streams, ponds, estuaries, oceans, or other waters unless in accordance with the requirements of a National Discharge Eliminations System (NDPES) permit and the permitting authority has been notified in writing prior to discharge. Do not discharge effluent containing this product to sewer systems without previously notifying the local sewage treatment plant authority. For guidance contact your State Water Board or Regional Office of the EPA.

PHYSICAL OR CHEMICAL HAZARDS

Ethylene oxide gas is extremely flammable. Do not use near flame, electrical sparks, or hot surfaces. Ground all equipment to prevent static sparks.

DIRECTIONS FOR USE

It is a violation of Federal law to use this product in a manner inconsistent with its labeling. It is the responsibility of the employer of any person engaged in the handling or application of this product to follow the requirements of 29 CFR 1910.1047. To be used only in an EOGas sterilizer for hospital, medical, and veterinary sterilization. Read package insert for complete sterilization instructions and additional precautions.

STORAGE AND DISPOSAL

Do not contaminate water, food, or feed by storage or disposal.

STORAGE: Store in a cool place and out of direct sunlight.


PESTICIDE DISPOSAL: Unwanted or expired EOGas cartridges should be returned to the manufacturer for disposal. Contact Andersen Products, Inc. (tel. no. 336-376-3000) for instructions. If unwanted or expired EOGas cartridges cannot be disposed of according to label instructions, contact your State Pesticide or Environmental Control Agency, or the Hazardous Waste Representative at the nearest EPA Regional Office for guidance.

CONTAINER DISPOSAL: Do not reuse empty box, used cartridge, or used liner bags. Wrap box, cartridge, and liner bags in paper and discard in the trash.

EOGas is a trademark of Andersen Products, Inc.
EPA Reg. No. 9417-NC-001 and 5,160,700.

P/N XXXX 012800

2/3/00, 11:28 AM

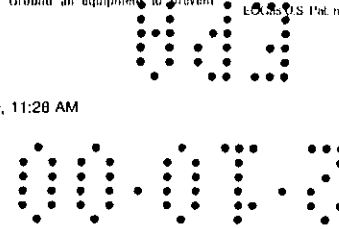

7/10



**AN
2011**



Use only in an Andersen EOGas Sterilizer and according to manufacturer's instructions. See outer box and package insert for precautions and sterilization instructions.

Active Ingredient: Ethylene Oxide	96%
Other Ingredient	4%
Total	100%

Net Contents: 0.39 av. oz. (11.00 g)

DANGER
KEEP OUT OF REACH OF CHILDREN
CONTENTS FLAMMABLE
DO NOT OPEN THIS PACKAGE
DISPOSE OF SPENT EOGas CARTRIDGE IN TRASH



LOT NO.

**USE
BEFORE**

Manufactured by: Andersen Sterilizers, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
Distributed by: Andersen Products, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
Distribué en Europe par: H. W. Andersen Products Ltd. • Clacton-On-Sea • Essex CO15 4XA U.K.
EPA Registration No. X00000-X • EPA Est. No. 9417-NC-001 • U.S. Patent No. 4,937,046 & 5,160,700


PN XXXX 012800

REMOVE TAPE AND TRIGGER GUARD, PRESS BUTTON TO ACTIVATE CARTRIDGE

ASSEMBLED
with 100% recycled
in EPA approved facility

MAY 10 2000

Under the Federal Insecticide,
Fungicide, and Rodenticide Act as
amended, for the pesticide,
registered under EPA Reg. No. *69340-4*

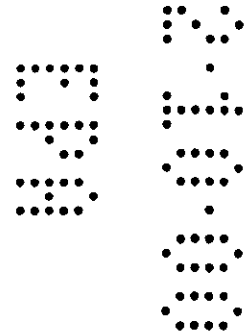

8/10


**AN
2014**

**KEEP OUT OF REACH OF CHILDREN
DANGER**

IMPORTANT: For Professional Use Only

Users must follow the requirements of the OSHA occupational exposure standard for Ethylene oxide (29 CFR 1910.1047)

Active ingredient: Ethylene oxide 96.0%
Inert ingredients: 4.0%
TOTAL 100.0%

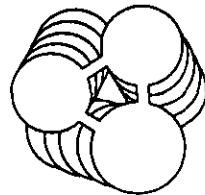

EOGas™
INDUSTRIAL REFILL PACK

ONE PACK

CONTENTS: 50 Cartridges. Each AN2014 cartridge contains 0.49 av. oz. (14.0 g) EOGas sterilant.

ANDERSEN STERILIZERS, INC.

Manufactured and Distributed by:
ANDERSEN STERILIZERS, INC.
3154 Caroline Drive • Haw River, NC 27258 USA
336-376-8622

Distributed by:
ANDERSEN PRODUCTS, INC.
3202 Caroline Drive • Haw River, NC 27258 USA
1-800-523-1276 • 336-376-3000
www.anpro.com • mailbox@anpro.com

European Distributor:
H. W. ANDERSEN PRODUCTS, LTD.
Clacton-On-Sea • Essex CO15 4XA UK

EPA Registration No. XXXXX-X
EPA Establishment No. 9417-NC-001



AN2014 Kit 012800

**ACCEPTED
with COMMENTS
in EPA Letter Dated:**

MAY 10 2000

Under the Federal Insecticide,
Fungicide, and Rodenticide Act as
amended, for the pesticide,
registered under EPA Reg. No. 69340-4

IF SWALLOWED: Drink at least two glasses of water. Do not induce vomiting. Do not give anything by mouth to an unconscious victim. Call a physician.

NOTE TO PHYSICIAN

Ethylene oxide is a gas. Skin exposure from contact with fabric, rubber or plastic containing residual ethylene oxide will commonly result in skin irritation with extensive blister formation. At high concentrations severe conjunctivitis can occur. Irritation of the respiratory tract may occur, but without acute lung edema. Symptoms of systemic intoxication are headache, nausea, vomiting, incoordination, and cardiac irregularities. Treatment is symptomatic.

ENVIRONMENTAL HAZARDS

Do not discharge effluent containing this product into lakes, streams, ponds, estuaries, oceans, or other waters unless in accordance with the requirements of a National Discharge Eliminations System (NDPES) permit and the permitting authority has been notified in writing prior to discharge. Do not discharge effluent containing this product to sewer systems without previously notifying the local sewage treatment plant authority. For guidance contact your State Water Board or Regional Office of the EPA.

PHYSICAL OR CHEMICAL HAZARDS

Ethylene oxide gas is extremely flammable. Do not use near flame, electrical sparks, or hot surfaces. Ground all equipment to prevent static sparks.

DIRECTIONS FOR USE

It is a violation of Federal law to use this product in a manner inconsistent with its labeling. It is the responsibility of the employer of any person engaged in the handling or application of this product to follow the requirements of 29 CFR 1910.1047. To be used only in an EOGas sterilizer for hospital, medical, and veterinary sterilization. Read package insert for complete sterilization instructions and additional precautions.

STORAGE AND DISPOSAL

Do not contaminate water, food, or feed by storage or disposal.
STORAGE: Store in a cool place and out of direct sunlight.


PESTICIDE DISPOSAL: Unwanted or expired EOGas cartridges should be returned to the manufacturer for disposal. Contact Andersen Products, Inc. (tel. no. 336-376-3000) for instructions. If unwanted or expired EOGas cartridges cannot be disposed of according to label instructions, contact your State Pesticide or Environmental Control Agency, or the Hazardous Waste Representative at the nearest EPA Regional Office for guidance.

CONTAINER DISPOSAL: Do not reuse empty box, used cartridge, or used liner bags. Wrap box, cartridge, and liner bags in paper and discard in the trash.

EOGas® is a trademark of Andersen Products, Inc.
EPA Reg. No. S. Par. no. 4,937,046 and 5,160,700.

P/N XXXX 012800

2/3/00, 11:28 AM

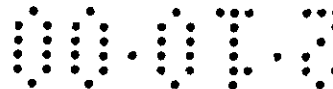

9/10



**AN
2014**



REMOVE TAPE AND TRIGGER GUARD, PRESS BUTTON TO ACTIVATE CARTRIDGE

Use only in an Andersen EOGas Sterilizer and according to manufacturer's instructions. See outer box and package insert for precautions and sterilization instructions.

Active Ingredient: Ethylene Oxide	96%
Other Ingredient	4%
Total.....	100%

Net Contents: 0.49 av. oz. (14.00 g)

LOT NO.

USE
BEFORE

DANGER
KEEP OUT OF REACH OF CHILDREN
CONTENTS FLAMMABLE
DO NOT OPEN THIS PACKAGE
DISPOSE OF SPENT EOGas CARTRIDGE IN TRASH



Manufactured by: Andersen Sterilizers, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
Distributed by: Andersen Products, Inc. • Health Science Park • Haw River, NC 27258 U.S.A.
Distribut en Europe par: H. W. Andersen Products Ltd. • Clacton-On-Sea • Essex CO15 4XA U.K.
EPA Registration No. XXXXX-X • EPA Est. No. 9417-NC-001 • U.S. Patent No. 4,937,048 & 5,160,700


PN XXXX 012800

REMOVE TAPE AND TRIGGER GUARD, PRESS BUTTON TO ACTIVATE CARTRIDGE

ACCEPTED
with **COMMENTS**
in EPA Letter Dated:

MAY 10 2000.

Under the Federal Insecticide,
Fungicide, and Rodenticide Act as
amended, for the pesticide,
registered under EPA Reg. No. 69340-4

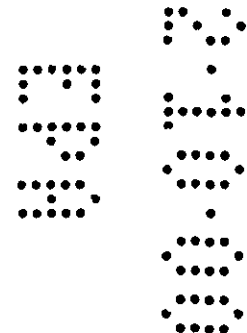

01/01