

70590-1

12-19-2002

1/2

{Front Side of Pouch}

WHEN OPENING: AVOID CONTACT WITH CLOTHING

←===== Tear Here =====>

Hype-Wipe®

DISINFECTING TOWEL WITH BLEACH

For One-Step Cleaning, Disinfecting Counters, Equipment and other Hard, Non-Porous Surfaces.

Contains at least a 1:10 dilution (0.55% by weight or 5500 ppm) of NaOCl (bleach) at expiration with proper storage.

KILLS:

Staphylococcus aureus • *Pseudomonas aeruginosa* • *Salmonella choleraesuis*


Active Ingredient:

Sodium Hypochlorite 0.94%

Other Ingredients 99.06%

Total 100.00%

Towel: Folded 6" x 12" nonwoven polypropylene
EPA Reg. No. 70590-1 EPA Est. 70590-IN-1
Net Wt. 12.2 grams (excludes towel)

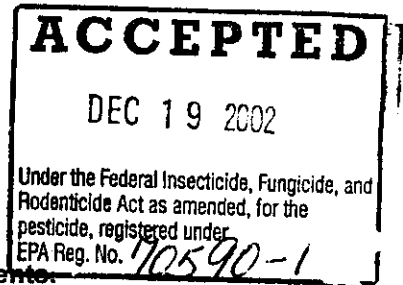

Keep Out of Reach of Children

CAUTION

See back panel for additional precautionary statements.

current technologies, inc.

P.O. Box 21

Crawfordsville, IN 47933 U.S.A.

Phone: 765-364-0490 Fax: 765-364-1607

Reorder # 9103

{Back Side of Pouch}

WHEN OPENING: AVOID CONTACT WITH CLOTHING

←===== Tear Here =====>

- Meets AOAC Hard Surface Carrier Disinfection Test
- Contains at least a 1:10 dilution of household bleach (5500 ppm) at expiration date (stamped on side [back] of pouch) if properly stored.
- Suitable for hard, non-porous surfaces including stainless steel, plastics, glass, ceramics, vinyl, porcelain.
- For use in clinical and industrial labs, clinics, hospitals, blood banks, physician and dental offices, wastewater facilities, emergency units, nursing stations.

This product is not to be used as a terminal sterilant/high level disinfectant on any surface or instrument that (1) is introduced directly into the human body, either into or in contact with the bloodstream or normally sterile areas of the body, or (2) contacts intact mucous membranes but which does not ordinarily penetrate the blood barrier, or otherwise enter normally sterile areas of the body. This product may be used to preclean or decontaminate critical or semi-critical medical devices prior to sterilization or high level disinfection.

2/2

DIRECTIONS FOR USE

It is a violation of Federal law to use this product in a manner inconsistent with its labeling.

1. Open pouch, remove towel. Use towel and liquid to wipe surface.
2. Allow solution to contact surface for 5 minutes before wiping or allow to air dry. Dispose towel in appropriate waste container.
3. **Personal Protection:** When handling surfaces soiled with blood/body fluids, wear disposable gloves, gowns and eye covering.

Storage: Store in a cool, dry area away from sunlight and heat to inhibit deterioration/reduced concentration of bleach solution. **Disposal:** Do not reuse empty pouch; wrap and put in trash.

PRECAUTIONARY STATEMENTS

Physical or Chemical Hazards

Do not use this product with ammonia, acids (such as vinegar), rust removers, toilet bowl cleaners. Will corrode aluminum.

Expiration Date: {MM/DD/YY}

[] Denotes alternate/optional language

{ } Denotes language that does not appear on the market label

